

YMC-Triart

Versatile hybrid silica based (U)HPLC columns

YMC
EUROPE GMBH
The Selectivity Company



Transfer

Scalable particles:
EASY
UHPLC ↔ HPLC


Flexible

YMC-Triart:
pH 1-12
Temperatures
up to 70°C

Universal


YMC-Triart
for acidic, basic and
neutral analytes

www.ymc.de

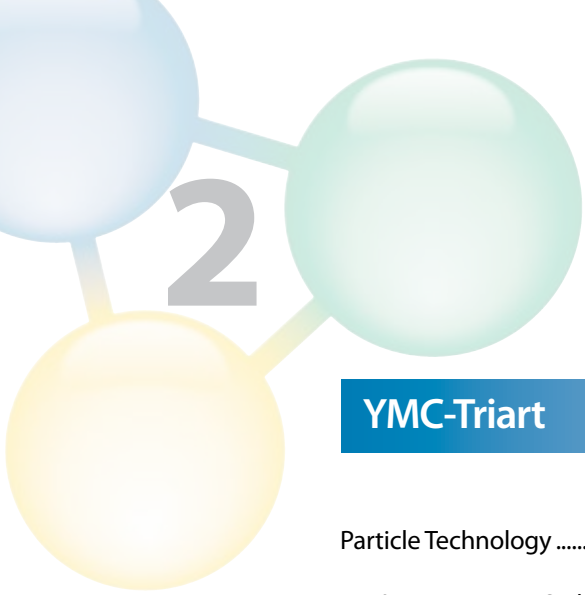


YMC-Triart

Particle Technology	page 4
pH & Temperature Stability	page 5
Method Transfer between HPLC ↔ UHPLC	page 6-7
UHPLC	page 8-19
LC/MS	page 20-22
Pharmaceuticals	page 23-28
Food	page 29-30
Peptides	page 31
HILIC	page 32-33
QC Data	page 34-39
Substance Index	page 40-41
Ordering Information	page 42-43

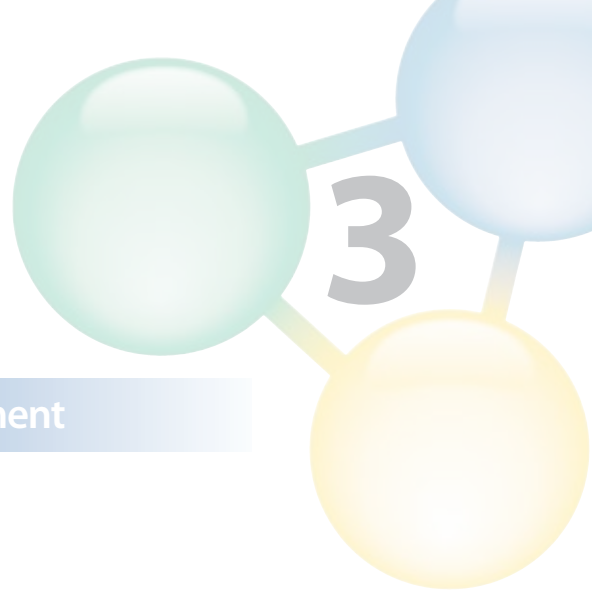

First choice column for method development

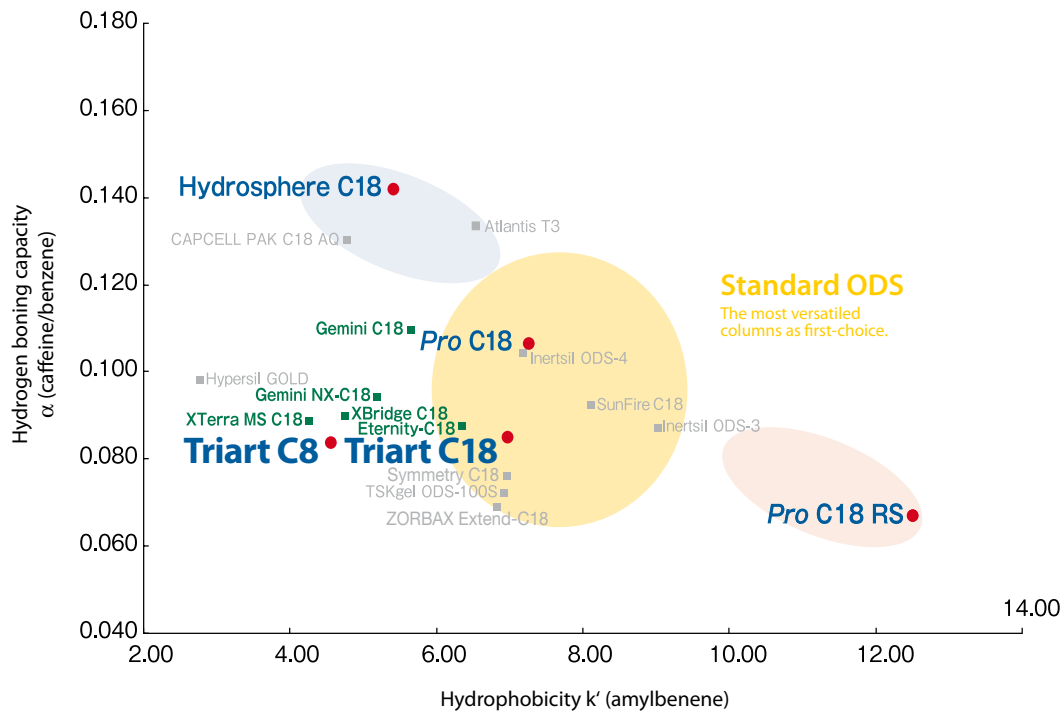

Conventional hybrid silica-based ODS columns tend to be less hydrophobic than silica-based columns. YMC-Triart C18 has a higher carbon load, giving it a hydrophobicity comparable to that of standard ODS columns, thereby making it a "versatile first-choice" column for method development.


Chromatographers always seek to push the limits of HPLC columns to greater extremes to allow them to perform day-to-day with ever-changing pH, buffers and temperature ranges. The column for the laboratory of today must be suitable for harsh pH conditions in combination with high temperature ranges without sacrificing selectivity. In addition narrow, symmetrical peak shapes are necessary in order to cope with rapid analysis of demanding samples. This has required manufacturers to seek more innovative ways to produce suitable stationary phases.

In order to meet these goals, YMC has developed a new particle technology. This is based on a multi-layered particle produced via a tightly controlled granulation technology which has been adapted from micro-reactor technology. The revolutionary production technique provides a multi-layer silica-organic hybrid stationary phase, which provides an outstandingly narrow pore size and particle size distribution. This in turn, results in low back pressures and high loadability.

Particle technology

YMC-Triart is a multi-layered material prepared using tightly controlled particle formation technology which has been adapted from micro-reactor technology. This recently developed production process results in exceptionally narrow particle and pore size distributions. With YMC-Triart, challenging pH and high temperature conditions are no longer a limitation to the day-to-day work in laboratories. Most importantly, due to its unique particle composition, a balanced hydrophobicity and silanol activity are achieved which makes YMC-Triart a "First Choice" column in method development.

YMC-Triart hybrid structure

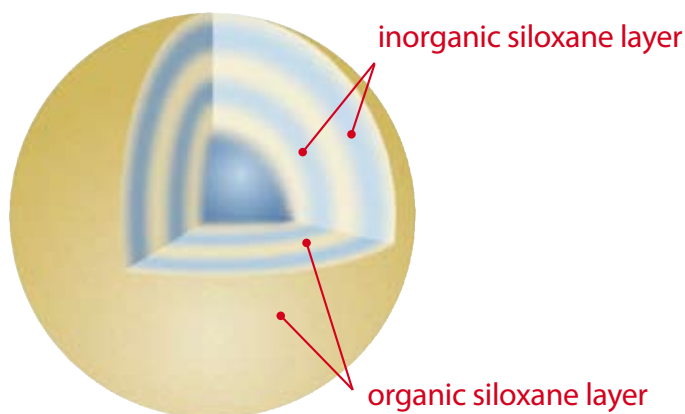


Specification

	YMC-Triart C18	YMC-Triart C8	Triart Phenyl	Triart PFP	YMC-Triart Diol-HILIC
Base	organic/inorganic silica				
Stationary phase	C18 (as USP L1)	C8 (as USP L7)	Phenyl (USP L11)	Pentafluorophenyl (USP L43)	Diol (USP L20)
Particle size	1.9, 3 and 5 μm				
Pore size	12 nm				
Bonding	polymeric type	polymeric type	monomeric type	monomeric type	monomeric type
End-capping	multi-stage hybrid groups	multi-stage hybrid groups	multi-stage hybrid groups	none	none
pH range	1 ~ 12	1 ~ 12	1 ~ 10	1 ~ 8	2 ~ 10
Temperature range	pH 1-7: 70 °C, pH 7-12: 50 °C	pH 1-7: 70 °C, pH 7-12: 50 °C	50 °C	50 °C	50 °C


Versatile wide pH stability

Phosphate buffer (pH 11.5, 40 °C)

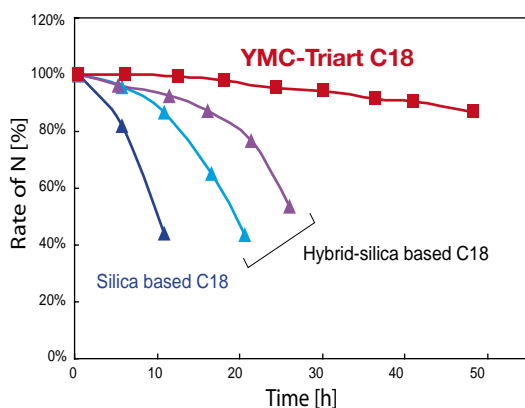

Column: 5 μ m, 150 x 4.6 mm ID
 Part-No.: TA12S05-1546WT
 Eluent: 50 mM K_2HPO_4 - K_3PO_4 (pH 11.5) / methanol (90/10)
 Flow rate: 1.0 ml/min
 Sample: benzyl alcohol


Triethylamine (pH 11.5, 40 °C)

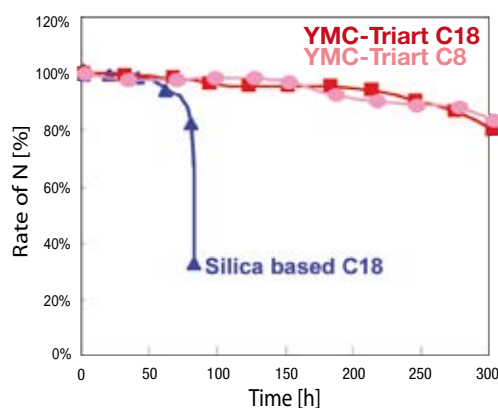

Column: 5 μ m, 150 x 4.6 mm ID
 Part-No.: TA12S05-1546WT
 Eluent: 50 mM triethylamine (pH 11.5) / methanol (90/10)
 Flow rate: 1.0 ml/min
 Sample: benzyl alcohol


Durability at high temperature

pH 6.9, 70 °C

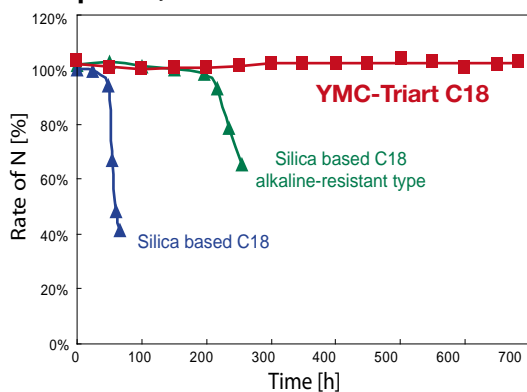

Column: 5 μ m, 50 x 2.0 mm ID
 Part-No.: TA12S05-0502WT
 Eluent: 20 mM KH_2PO_4 - K_2HPO_4 (pH 6.9) / acetonitrile (90/10)
 Flow rate: 0.2 ml/min
 Temperature: 70 °C
 Sample: phenol

pH 1, 70 °C

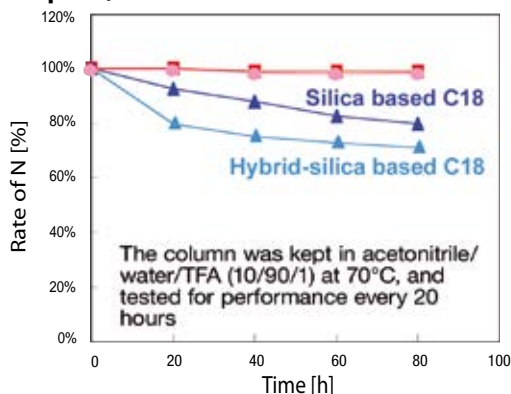

Column: 5 μ m, 50 x 2.0 mm ID
 Part-No.: TA12S05-0502WT
 Eluent: acetonitrile / water (60/40)
 Flow rate: 0.2 ml/min
 Temperature: 70 °C
 Sample: butyl benzoate

YMC-Triart phases show great chemical stability due to the newly developed hybrid-silica. Even under high pH or high temperature conditions, the lifetime of YMC-Triart phases is more than 10x greater than conventional reversed phase columns.

Application Data by courtesy YMC Co., Ltd.


6

Transfer
HPLC
↕
UHPLC

Secure your method transfer!

Differences in selectivity, retention time, and also peak shapes between different particle sizes of commercially available C18 phases in the same brand (or an alternative as recommended by its manufacture) have been observed.


YMC-Triart C18

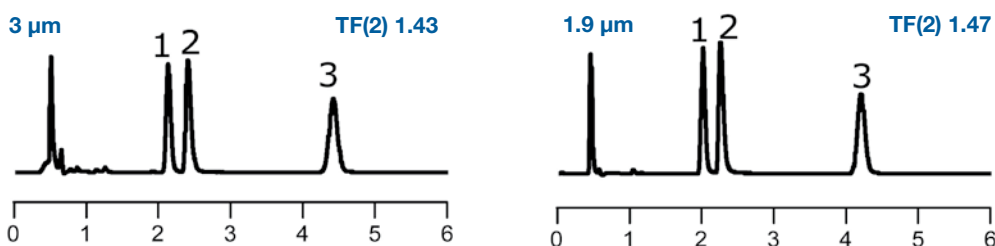


YMC has addressed this issue of method transfer. YMC-Triart columns show identical selectivity and excellent peak shapes for basic compounds for all 3.0 µm to 1.9 µm particle sizes. It allows predictable scale up from UHPLC to conventional HPLC and even to semi-preparative LC, and vice versa.

Case Studies*

X-Bridge BEH C18 and Acquity UPLC BEH C18

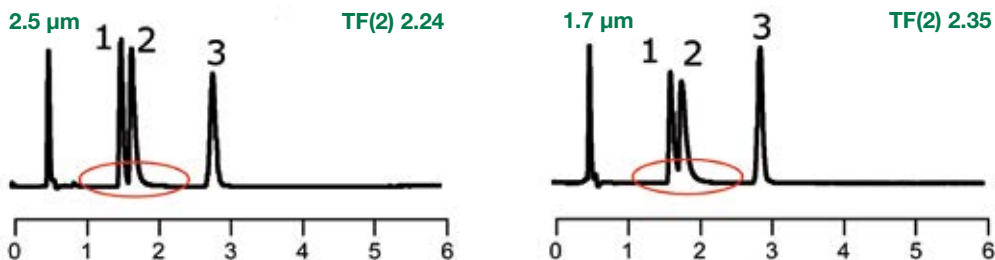

Kinetex™ C18

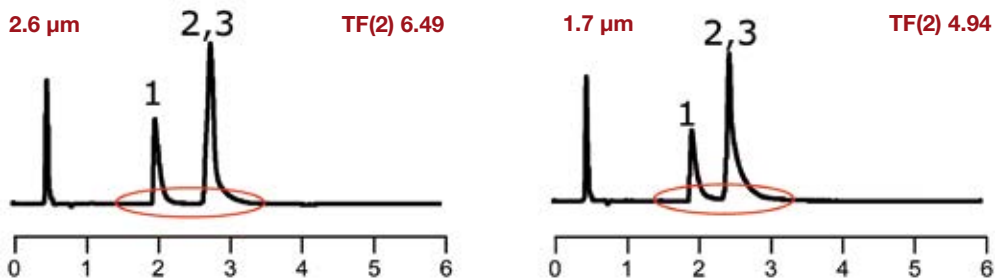


Kinetex™ C18 columns show significant peak tailing™ and have limited scalability due to lack of larger particle sizes.

Column: 50 x 2.0 mm ID or 2.1 mm ID
Eluent: 20 mM KH_2PO_4 - K_2HPO_4 (pH 6.9) / acetonitrile (65/35)
Temperature: 40 °C
Flow rate: 0.2 ml/min
Detection: UV at 235 nm

1. Chlorpheniramine (basic)
2. Dextromethorphan (basic)
3. Propyl paraben (internal standard)


* These observations might not be representative for all applications but have been reported in some cases.

Evaluation of method transfer performance!

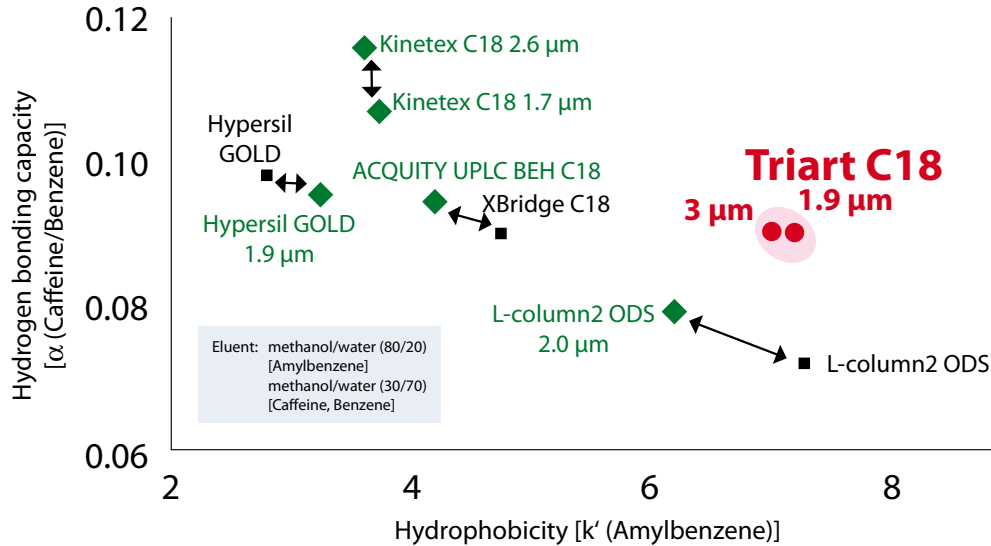

With the introduction of UHPLC, sub-2- μ m particles became necessary. Therefore smaller particles have been added to existing column lines. Consequently, sub-2- μ m particles may exhibit differences in chromatographic performance. By introducing YMC-Triart, YMC provides matching chromatographic behaviour for **all** particles sizes!

Method transfer between HPLC ↔ UHPLC

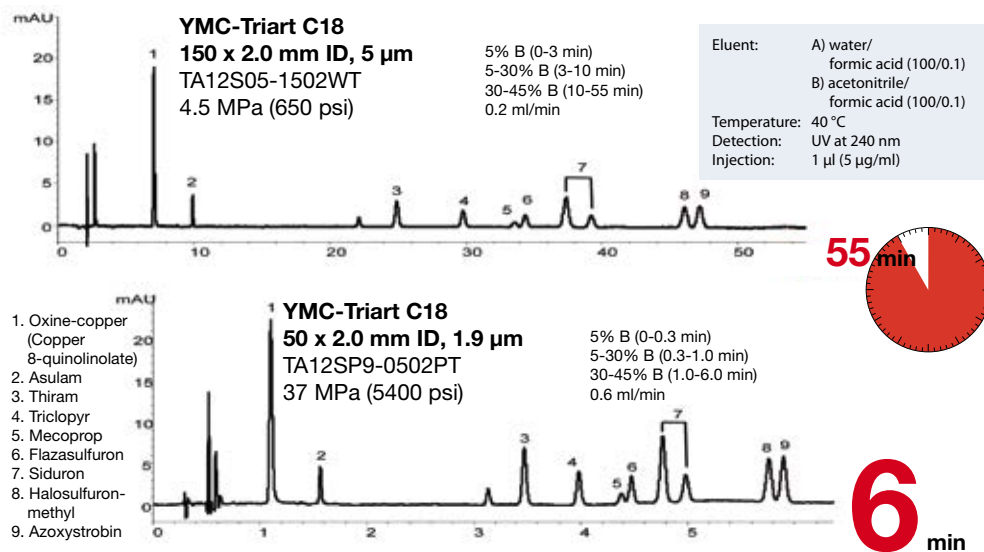


Application Data by courtesy YMC Co., Ltd.

8


UHPLC

Higher resolution and good loadability




YMC-Triart C18

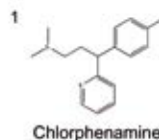

Ascentis Express C18

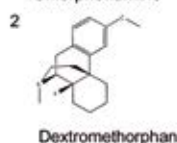


Kinetex C18

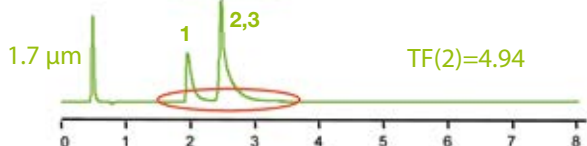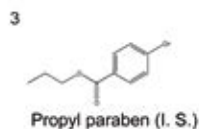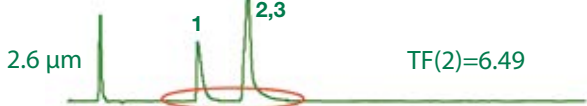

Column: 50 x 2.0 or 2.1 mm ID
 Eluent: 20 mM KH₂PO₄-K₂HPO₄ (pH 6.9)/ acetonitrile (65/35)
 Flow rate: 0.2 ml/min
 Detection: UV at 235 nm
 Temperature: 40 °C




Lower HETP means higher resolution!

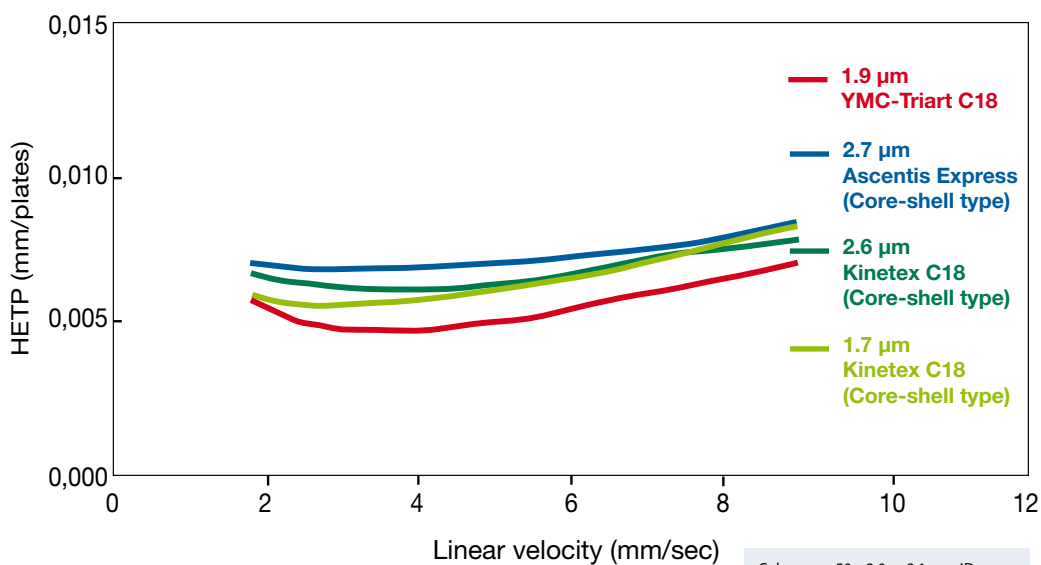


Column: 50 x 2.0 or 2.1 mm ID
 Eluent: acetonitrile / water (60/40)
 Detection: UV at 254 nm
 Sample: Butylbenzoate

Application Data by courtesy YMC Co., Ltd.


Determination of Artificial Sweeteners with LC-MS/MS

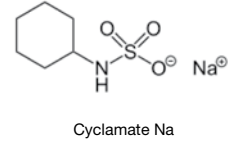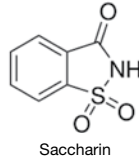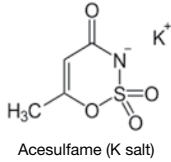


→ Non biological markers of wastewater entries in ground and surface water

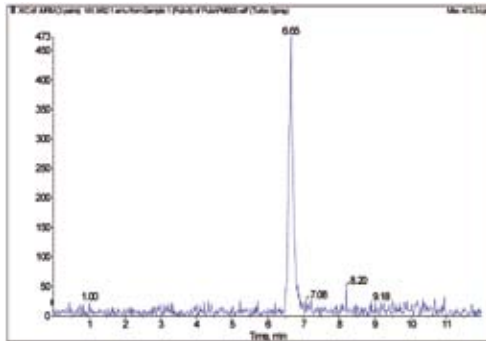

Extracted Ion Chromatogram (XIC) of Acesulfame K, 0.1 µg/L

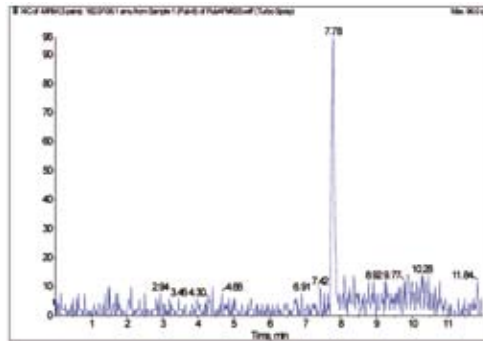

Extracted Ion Chromatogram (XIC) of Saccharin, 0.1 µg/L

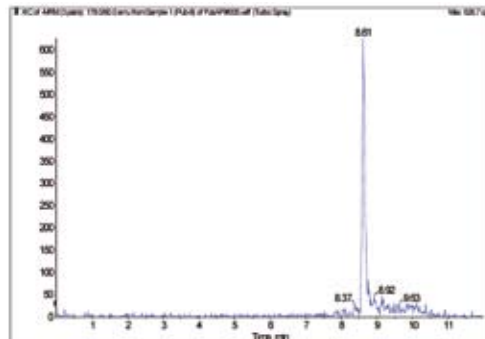

Extracted Ion Chromatogram (XIC) of Cyclamate Na, 0.1 µg/L


Column: YMC-Triart C18, 12 nm, 1.9 µm, 100 x 3.0 mm ID
 Part-No.: TA12SP9-1003PT
 LC-System: Agilent 1100 HPLC system and CTC Analytics
 HTC-Pal Autosampler
 MS/MS System: Applied Biosystems MDS Sciex API 4000,
 ESI negative
 Temperature: 35°C
 Flow: 0.3 ml/min
 Injection: 40 µL, direct injection
 Eluent: A: H₂O (containing 10 mmol NH₄ formate)
 B: MeOH (containing 10 mmol NH₄ formate)
 Gradient:

Time	0	6.0	6.1	12.0
% B	2	75	2	2

by courtesy of: Thomas Class, Sandro Jooß
 PTRL Europe, Helmholtzstraße 22, Science Park I, D-89081 Ulm


10 UHPLC

Angiotensin I, II and III

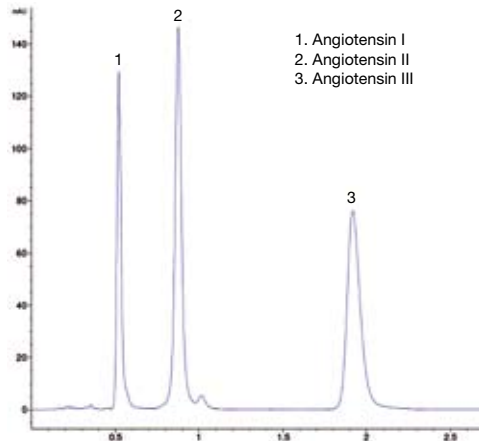


Column: YMC-Triart C18, 1.9 μ m (50 x 2.0 mm ID)
Part-No.: TA12SP9-0502PT
Eluent: 20 mM KH_2PO_4 + K_2HPO_4 (pH 7.9) / acetonitrile (22/78)
Flow rate: 0.7 ml/min
Detection: UV at 220 nm
Pressure: 720 bar
Injection: 0.5 μ l
Temperature: 40 $^\circ\text{C}$

Insulin

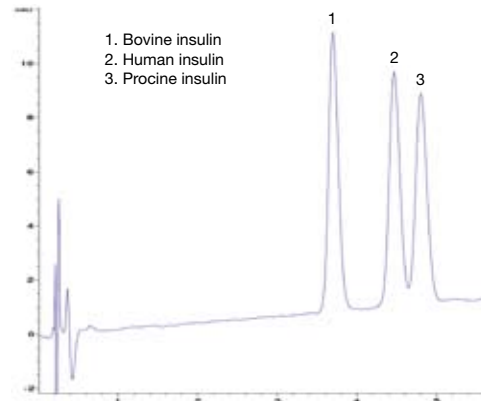


Column: YMC-Triart C18, 1.9 μ m (50 x 2.0 mm ID)
Part-No.: TA12SP9-0502PT
Eluent: A) H_2O + 0.1% TFA
B) acetonitrile + 0.1% TFA
Gradient: 30% B (0 min); 30-32% B (0-5 min); 32% B (55 min)
Flow rate: 0.6 ml/min
Detection: UV at 220 nm
Pressure: 611 bar
Injection: 0.5 μ l
Temperature: 30 $^\circ\text{C}$

Macrolide antibiotics

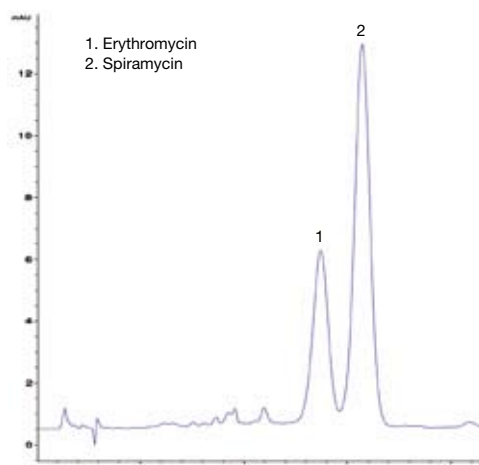


Column: YMC-Triart C18, 1.9 μ m (50 x 2.0 mm ID)
Part-No.: TA12SP9-0502PT
Eluent: A) 20 mM K_2HPO_4 + 20 mM KH_2PO_4 (pH 7.9)
B) acetonitrile
Gradient: 60% B (0.5 min); 60-70% B (0.5-1.5 min); 70% B (3.5 min)
Flow rate: 0.45 ml/min
Detection: UV at 210 nm
Pressure: 520 bar
Injection: 1 μ l
Temperature: 50 $^\circ\text{C}$

Sulpha drugs

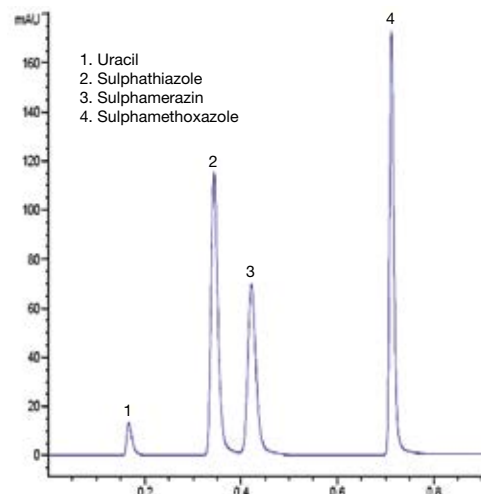


Column: YMC-Triart C18, 1.9 μ m (50 x 2.0 mm ID)
Part-No.: TA12SP9-0502PT
Eluent: H_2O + formic acid (pH 2.5) / acetonitrile (75/25)
Flow rate: 0.75 ml/min
Detection: UV at 280 nm
Pressure: 740 bar
Injection: 0.5 μ l
Temperature: 50 $^\circ\text{C}$

Tetracycline antibiotics

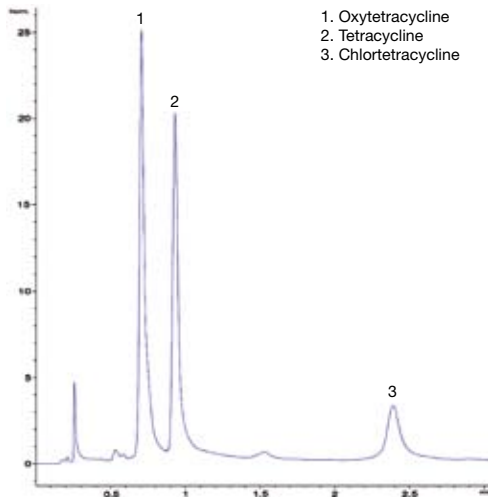


Column: YMC-Triart C18, 1.9 μm (50 x 2.0 mm ID)
Part-No.: TA12SP9-0502PT
Eluent: 5 mM $\text{NH}_4\text{CH}_2\text{COOH}$ / acetonitrile (87/13)
Flow rate: 0.65 ml/min
Detection: UV at 280 nm
Pressure: 662 bar
Injection: 1 μl
Temperature: 40 $^\circ\text{C}$

Betablockers

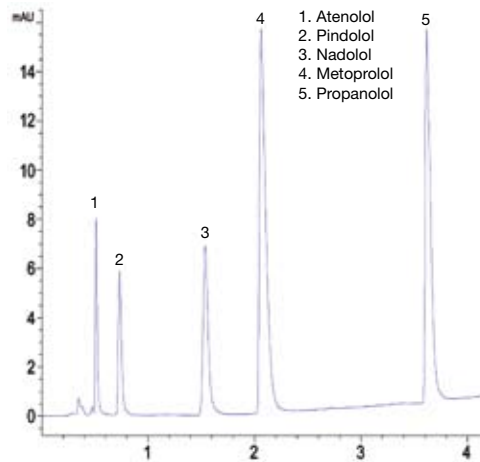

Column: YMC-Triart C18, 1.9 μm (50 x 2.0 mm ID)
Part-No.: TA12SP9-0502PT
Eluent: A) 20 mM $\text{NH}_4\text{CH}_2\text{COOH}$ + ammonia (pH 9.0)
B) acetonitrile
Gradient: 25% B (1.0 min); 75% B (1-6 min)
Flow rate: 0.35 ml/min
Detection: UV at 254 nm
Pressure: 450 bar
Injection: 1 μl
Temperature: 40 $^\circ\text{C}$

Fast LC for conventional HPLC

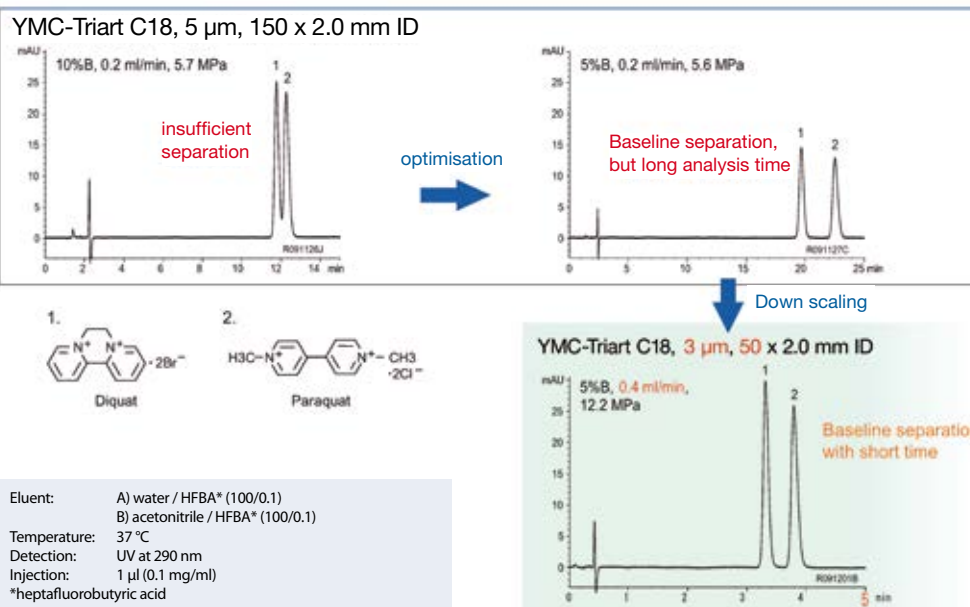


Application Data by courtesy YMC Co., Ltd.

12

Triart
Phenyl



Fexofenadine hydrochloride (Japanese Pharmacopoeia)

A) Standard solution *1 (0.001 mg/ml fexofenadine hydrochloride)

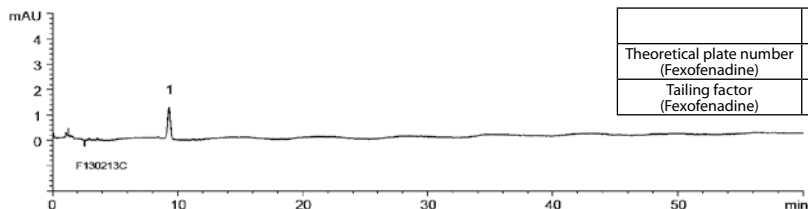

	System suitability requirement	result
Theoretical plate number (Fexofenadine)	≥ 8000	10100
Tailing factor (Fexofenadine)	≤ 2.0	1.00



B) Sample solution *1 (1 mg/ml fexofenadine hydrochloride)

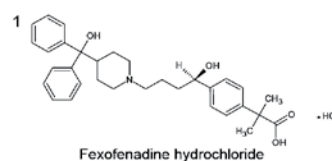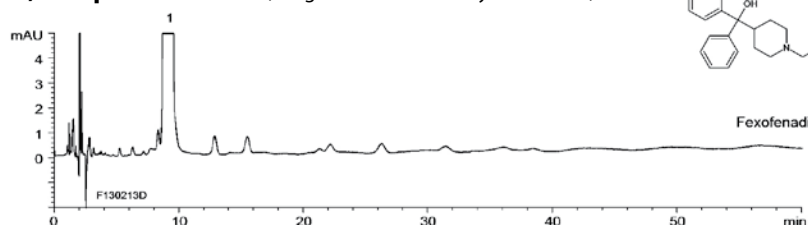

Column: YMC-Triart Phenyl, 5 μm, 12 nm (250 x 4.6 mm ID)
 Part-No.: TPH12505-2546WT
 Eluent: acetonitrile / buffer *2 / triethylamine (350/650/3)
 *2 Dissolve 7.51 g of NaH₂PO₄·2H₂O and 0.96 g of NaClO₄·H₂O in 1000 ml water, adjust pH 2.0 with H₃PO₄
 Flow rate: 2.0 ml/min (adjust the flow rate so that the retention time of fexofenadine is about 9 min)
 Detection: UV at 220 nm
 Injection: 20 μl
 Temperature: 25 °C
 (The Japanese Pharmacopoeia 16th; related substances)

*1 All standard and sample solutions were prepared from fexofenadine hydrochloride supplied as a reagent for laboratory use.

Standard solution *1 (0.06 mg/ml fexofenadine hydrochloride)

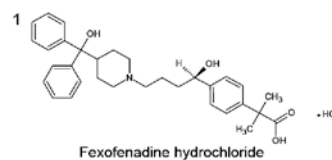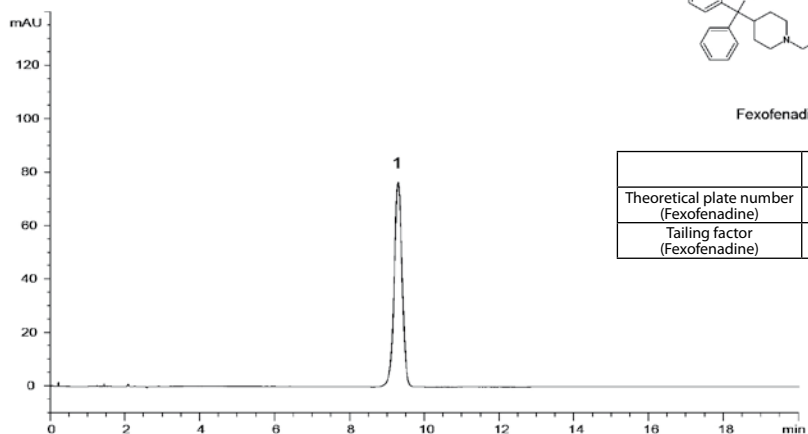


	System suitability requirement	result
Theoretical plate number (Fexofenadine)	≥ 8000	9500
Tailing factor (Fexofenadine)	≤ 2.0	0.98

Column: YMC-Triart Phenyl, 5 μm, 12 nm (250 x 4.6 mm ID)
 Part-No.: TPH12505-2546WT
 Eluent: acetonitrile / buffer *2 / triethylamine (350/650/3)
 *2 Dissolve 7.51 g of NaH₂PO₄·2H₂O and 0.96 g of NaClO₄·H₂O in 1000 ml water, adjust pH 2.0 with H₃PO₄
 Flow rate: 2.0 ml/min (adjust the flow rate so that the retention time of fexofenadine is about 9 min)
 Detection: UV at 220 nm
 Injection: 20 μl
 Temperature: 25 °C
 (The Japanese Pharmacopoeia 16th; assay)

*1 Standard solutions was prepared from fexofenadine hydrochloride supplied as a reagent for laboratory use.



Fexofenadine hydrochloride (US Pharmacopoeia)

A) Assay preparation *1 (assay), Reference solution *1 (related compounds) (0.06 mg/ml fexofenadine hydrochloride)

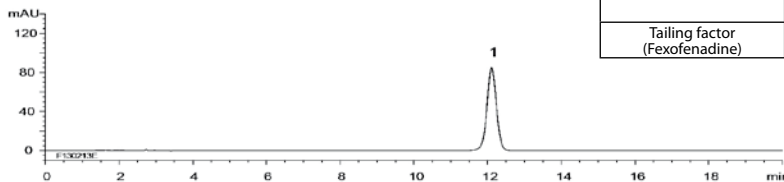

	System suitability requirement (assay)	result
Tailing factor (Fexofenadine)	≤ 2.0	1.00


B) Test solution *1 (related compounds) (1 mg/ml fexofenadine hydrochloride)

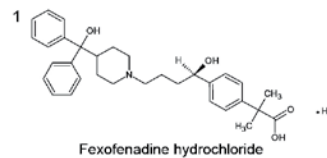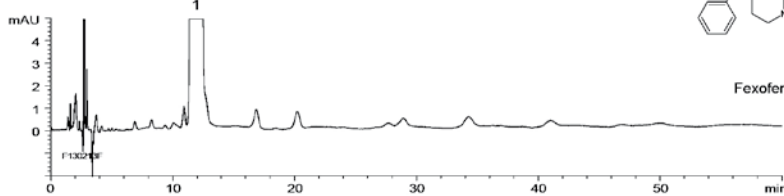

Column: YMC-Triart Phenyl, 5 µm, 12 nm (250 x 4.6 mm ID)
 Part-No.: TPH12S05-2546WT
 Eluent: acetonitrile / buffer *2 / triethylamine (350/650/3)
 *2 Dissolve 7.51 g of NaH₂PO₄·2H₂O and 0.96 g of NaClO₄·H₂O in 1000 ml water, adjust pH 2.0 with H₃PO₄
 Flow rate: 1.5 ml/min
 Detection: UV at 220 nm
 Injection: 20 µl
 Temperature: 25 °C
 (The United States Pharmacopoeia 36th; assay, related compounds)



*1 All standard and sample solutions were prepared from fexofenadine hydrochloride supplied as a reagent for laboratory use.

A) Standard solution *1 (0.015 mg/ml fexofenadine hydrochloride)

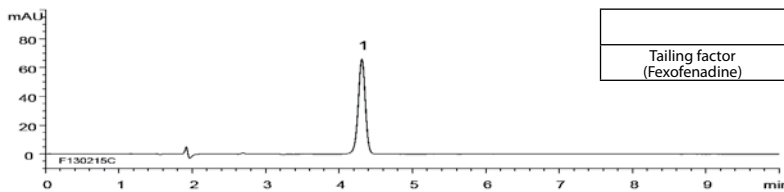

	System suitability requirement (assay)	result
Tailing factor (Fexofenadine)	≤ 2.0	0.95

B) Sample solution *2 (0.018 mg/ml fexofenadine hydrochloride)

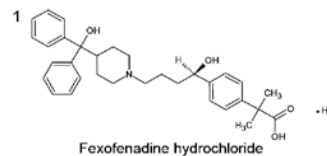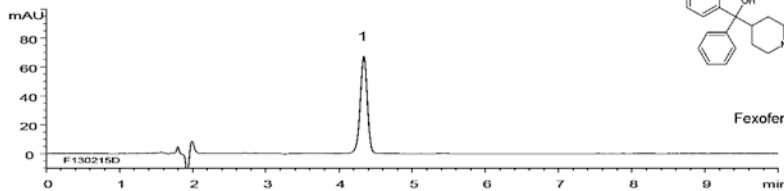


Column: YMC-Triart Phenyl, 5 µm, 12 nm (250 x 4.6 mm ID)
 Part-No.: TPH12S05-2546WT
 Eluent: acetonitrile / buffer *3 (9/16)
 *3 Add 15 ml of acetonitrile/triethylamine (1/1) to 1000 ml of acetic acid/water (17/9983), adjust pH 5.25 with H₃PO₄
 Flow rate: 1.5 ml/min
 Detection: UV at 220 nm
 Injection: 20 µl
 Temperature: 35 °C
 (The United States Pharmacopoeia 36th; assay)

*1 Standard solution was prepared from fexofenadine hydrochloride supplied as a reagent for laboratory use.



*2 Sample solution was prepared from fexofenadine hydrochloride tablets.

Analysis of coal tar dye (Brilliant blue FCF) and its impurities

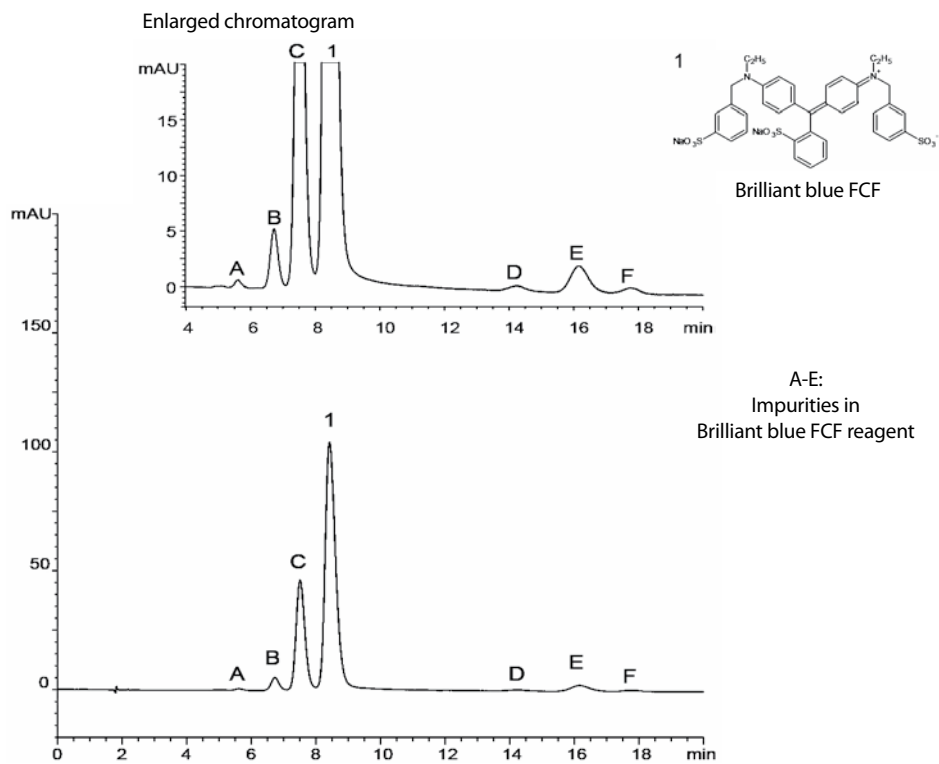


Column: YMC-Triart Phenyl, 5 μ m, 12 nm (150 x 3.0 mm ID)
 Part-No.: TPH12505-1503WT
 Eluent: 0.1% phosphoric acid / methanol (55/45)
 Flow rate: 0.425 ml/min
 Detection: UV at 630 nm
 Injection: 2.0 μ l
 Temperature: 40 $^{\circ}$ C


Catecholamines, serotonin, and their precursors and metabolites

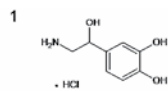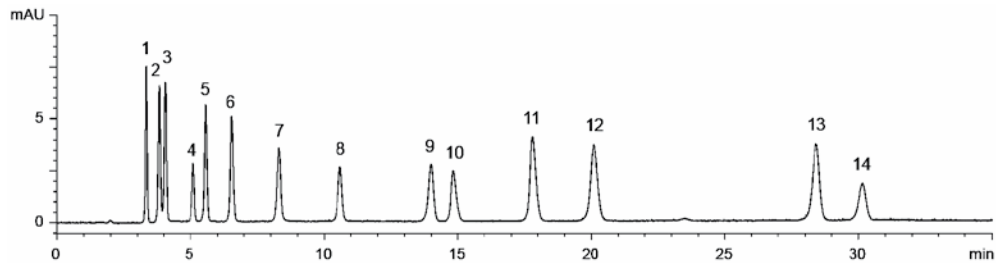


Noradrenaline hydrochloride (NA)
(Norepinephrine hydrochloride)

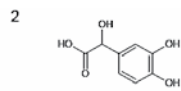


3,4-Dihydroxymandelic acid
(DOMA)

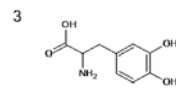


3,4-Dihydroxyphenylalanine
(DOPA)

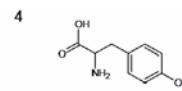


Tyrosine (Tyr)

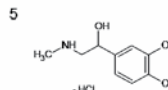


Adrenaline hydrochloride (A)
(Epinephrine hydrochloride)

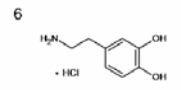


Dopamine hydrochloride (DA)

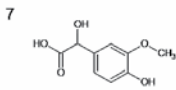


Vanillylmandelic acid (VMA)

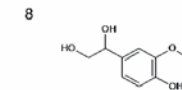


3-Methoxy-4-hydroxyphenylglycol
(MHPG)

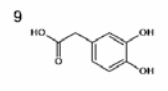


3,4-Dihydroxyphenylacetic acid
(DOPAC)

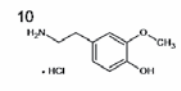


3-Methoxytyramine hydrochloride
(3MT)

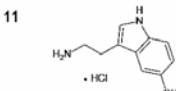


Serotonin hydrochloride
(5-Hydroxytryptamine hydrochloride, 5HT)

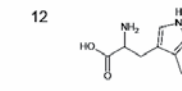

Tryptophan (Trp)

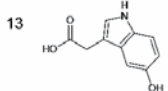

5-Hydroxyindoleacetic acid
(5HIAA)

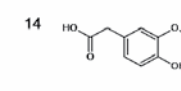


Homovanillic acid (HVA)

Column: YMC-Triart PFP, 3 μ m, 12 nm (150 x 3.0 mm ID)
 Part-No.: TPF12503-1503WT
 Eluent: A) 10 mM formic acid
 B) methanol containing 10 mM formic acid
 0-20% B (0-30 min), 20% B (30-35 min)
 Flow rate: 0.425 ml/min
 Detection: UV at 280 nm
 Injection: 4 μ l (5 μ g/ml)
 Temperature: 25 $^{\circ}$ C


16 UHPLC

Antidepressants

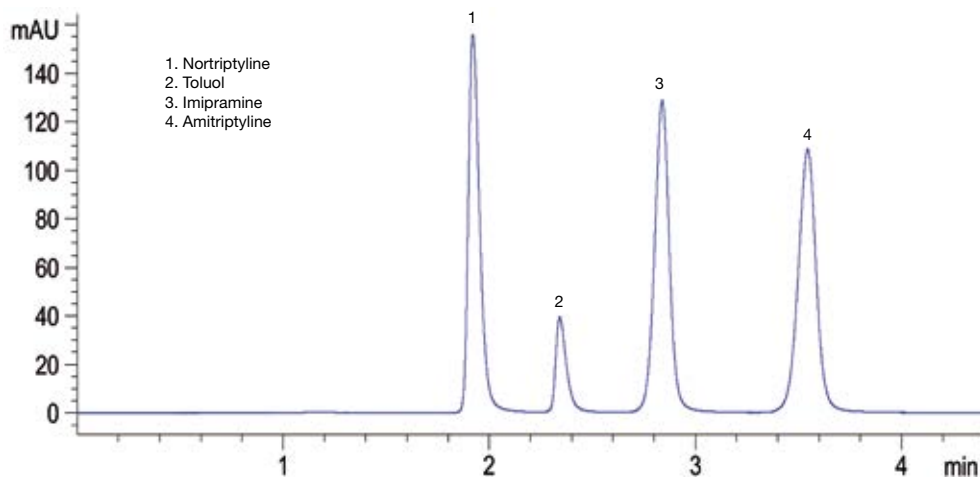


Column: YMC-Triart Phenyl, 1.9 μm (100 x 2.0 mm ID)
Part-No.: TPH12SP9-1002PT
Eluent: methanol / 25 mM KH_2PO_4 (pH 6.0) (65/35)
Flow rate: 0.4 ml/min
Detection: UV at 254 nm
Injection: 2 μl
Temperature: 25 $^\circ\text{C}$

Resorcinol

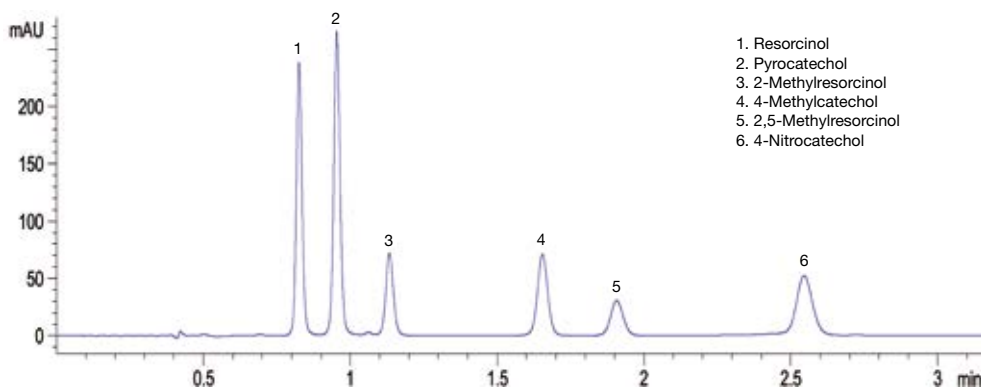


Column: YMC-Triart PFP, 1.9 μm (100 x 2.0 mm ID)
Part-No.: TPF12SP9-1002PT
Eluent: 0.1% formic acid / acetonitrile + 0.1% formic acid (85/15)
Flow rate: 0.8 ml/min
Detection: UV at 270 nm
Injection: 0.5 μl
Temperature: 25 $^\circ\text{C}$

Stevia leaves

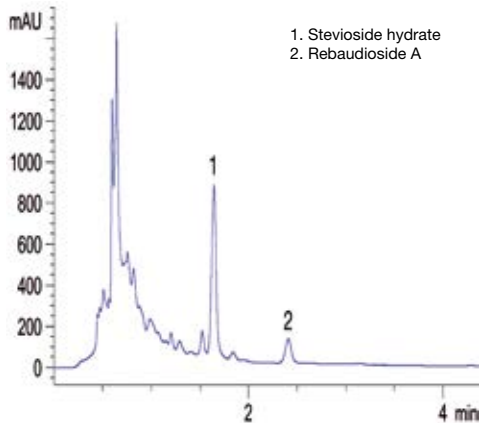

Column: YMC-Triart Diol-HILIC, 1.9 μm (100 x 3.0 mm ID)
Part-No.: TDH12SP9-1003PT
Eluent: acetonitrile / water (85/15)
Flow rate: 1 ml/min
Detection: UV at 200 nm
Injection: 2 μl
Temperature: 30 $^{\circ}\text{C}$

Nasal Spray

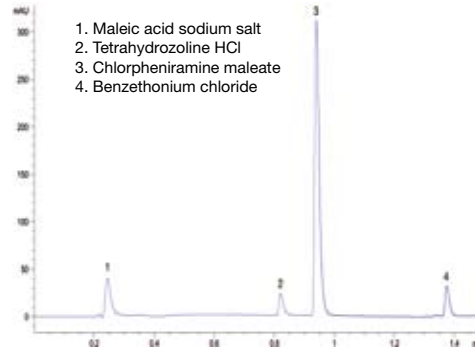


Column: YMC-Triart C18, 1.9 μm , 12 nm (50 x 2.0 mm ID)
Part-No.: TA12SP9-0502PT
Eluent: A) trifluoro acetic acid 0,05% / B) methanol (50/50)
Gradient:

min	A	B
0	80	20
0.5	10	90
1.2	0	100

Flow rate: 0.6 ml/min
Detection: UV at 260 nm
Injection: 0.2 μl
Temperature: 40 $^{\circ}\text{C}$

Eye drop formulation

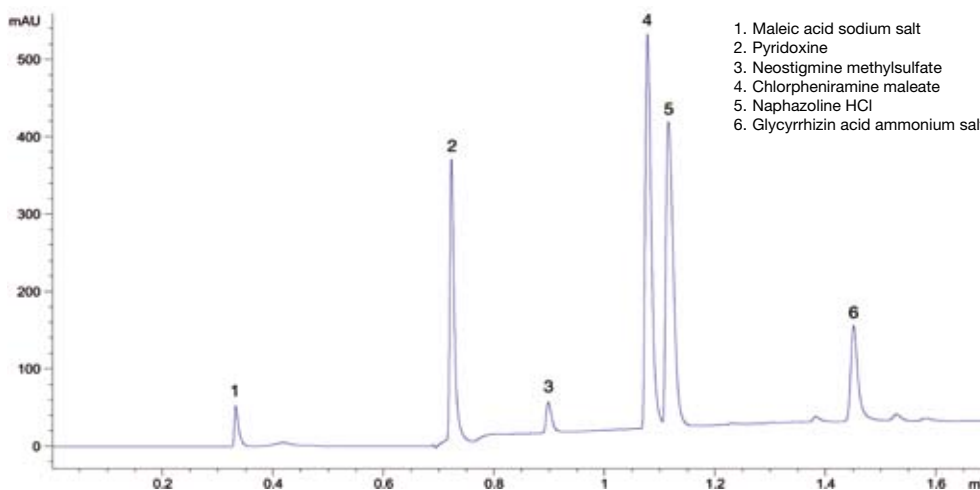


Column: YMC-Triart C18, 1.9 μm , 12 nm (50 x 2.0 mm ID)
Part-No.: TA12SP9-0502PT
Eluent: A) trifluoroacetic acid (0.05%)
B) acetonitrile
Flow rate: 0.6 ml/min
Detection: UV at 265 nm
Injection: 0.5 μl
Temperature: 40 $^{\circ}\text{C}$

Gradient:

min	A	B
0	100	0
1	50	50
1.5	50	50
1.7	10	90


Application Data by courtesy YMC Co., Ltd.

Guaiacol and impurities

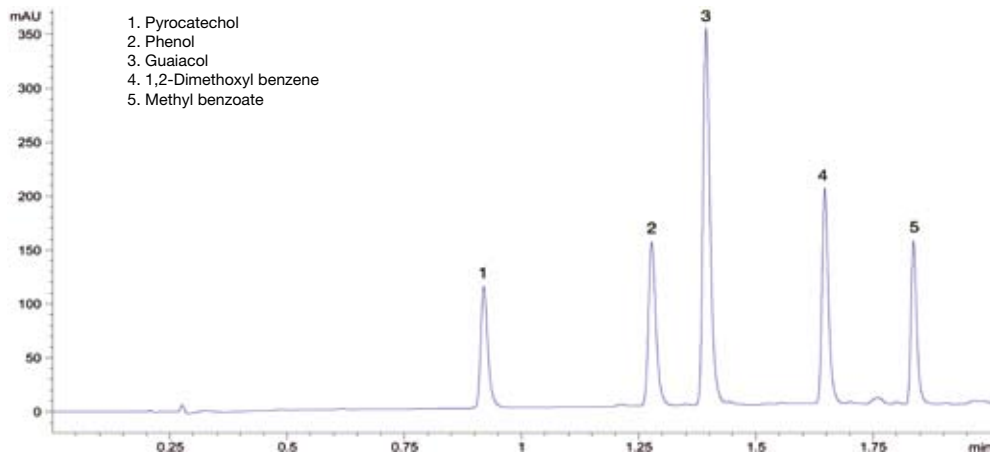

Column: YMC-Triart C18, 1.9 μ m, 12 nm (50 x 2.0 mm ID)
Part-No.: TA12SP9-0502PT
Eluent: water / acetonitrile (50/50)
Flow rate: 0.7 ml/min
Detection: UV at 254 nm
Injection: 0.5 μ l
Temperature: 40 $^{\circ}$ C


Paracetamol

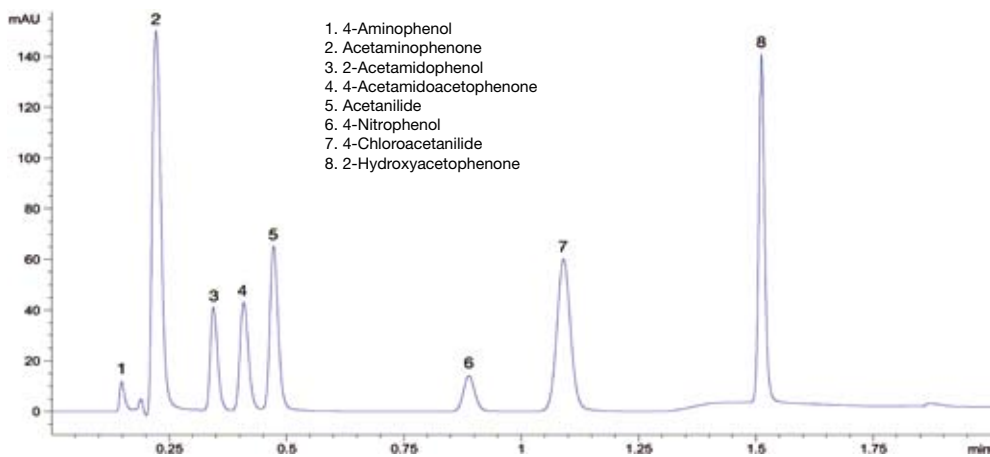

Column: YMC-Triart C18, 1.9 μ m, 12 nm (50 x 2.0 mm ID)
Part-No.: TA12SP9-0502PT
Eluent: A) formic acid / B) acetonitrile
Flow rate: 0.7 ml/min
Detection: UV at 254 nm
Injection: 0.5 μ l
Temperature: 40 $^{\circ}$ C

Gradient:	min	A	B
	0	70	30
	1	70	30
	1.5	20	80
	2	20	80

7 Analgesics

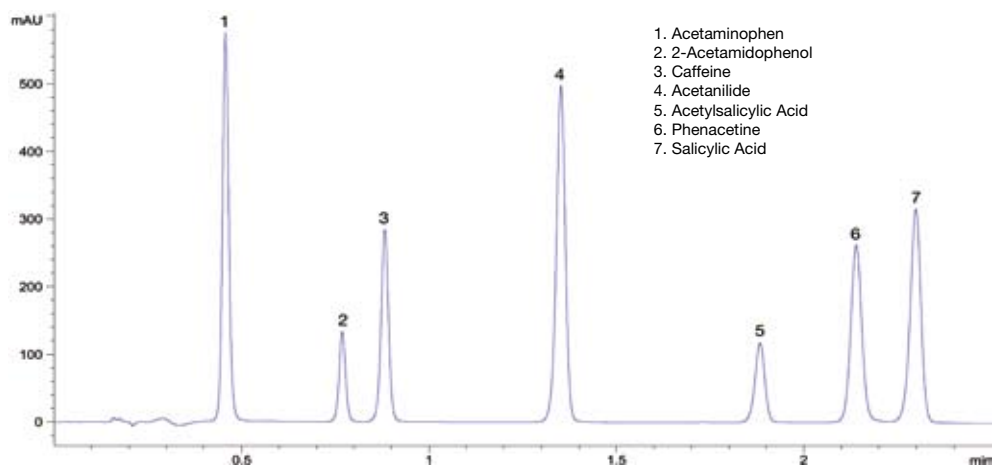


1. Acetaminophen
2. 2-Acetamidophenol
3. Caffeine
4. Acetanilide
5. Acetylsalicylic Acid
6. Phenacetine
7. Salicylic Acid

Column: YMC-Triart C18, 1.9 μ m, 12 nm (50 x 2.0 mm ID)
Part-No.: TA12SP9-0502PT
Eluent: formic acid / acetonitrile (50/50)
Flow rate: 0.8 ml/min
Pressure: 610 bar
Detection: UV at 240 nm
Injection: 1 μ l
Temperature: 40 $^{\circ}$ C

Application Data by courtesy YMC Co., Ltd.

20 LC/MS


LC/MS compatibility

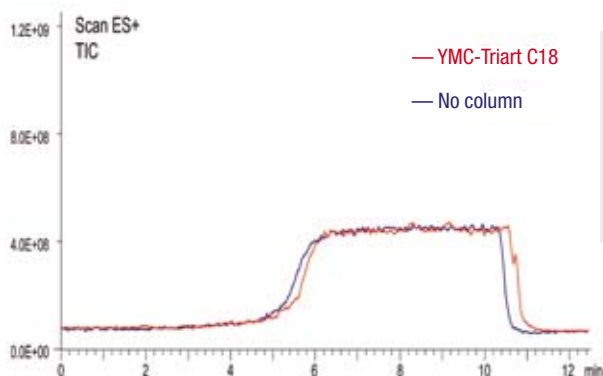

Column: 5 μ m, 50 x 2.0 mm ID
Part-No.: TA12S05-0502WT
Eluent: A) water / formic acid (100/0.1)
 B) acetonitrile / formic acid (100/0.1)
 5% B (0-1 min), 5-100% B (1-5 min),
 100% B (5-10 min), 100-5% B (10-10.1 min),
 5% B (10.1-12.5 min)
Flow rate: 0.4 ml/min
Temperature: 40 °C
Detection: ESI positive, TIC (Mass Range: 50-1000)


Column bleeding, caused by the fragments of stationary phase, is the main reason for background noise and restrictions on detection limits. No bleed is observed in the test of total ion current (TIC) measured by LC/MS with blank or with YMC-Triart C18. So in terms of the signal/noise ratio (S/N ratio), YMC-Triart C18 can be expected to not only reduce the background noise but to also increase the sensitivity of the analysis.

LC/MS analysis of benzodiazepine derivatives

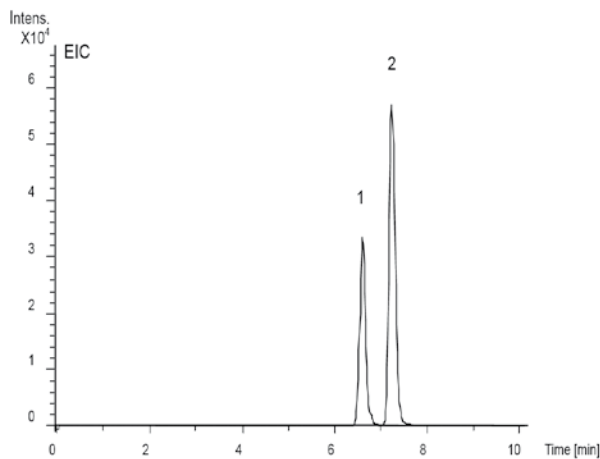

Peak 1: 100ng_mL_RB4_01_2291.d: EIC 343.051 \pm 0.01+All MS
 Peak 2: 100ng_mL_RB4_01_2291.d: EIC 343.078 \pm 0.01+All MS






Courtesy of J. Watanabe, Bruker Daltonics K. K.

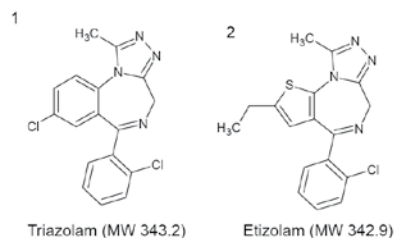



Column: YMC-Triart C18 (5 μ m, 12 nm)
 50 x 2.0 mm ID
Part-No.: TA12S05-0502WT
Eluent: A) 10 mM formic acid
 B) acetonitrile
Gradient: 25-50% B (0-10 min)
Flow rate: 0.2 ml/min
Temperature: 40 °C
Detection: Bruker Daltonics micrOTOF, ESI,
 positive mode
Injection: 5 μ l (100 ng/ml)

LC/MS analysis of tetracycline antibiotics

A) UV

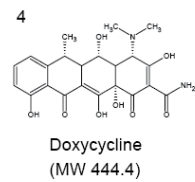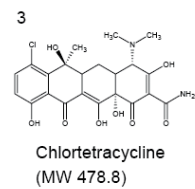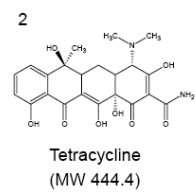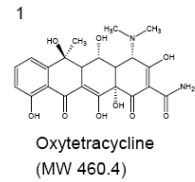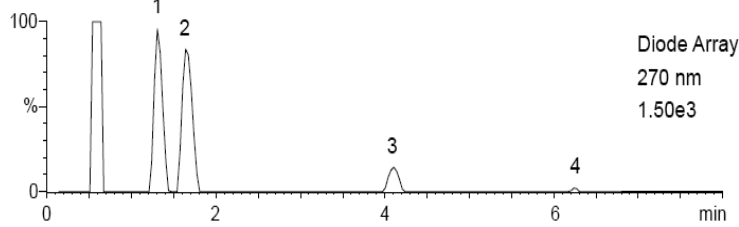

B) MS

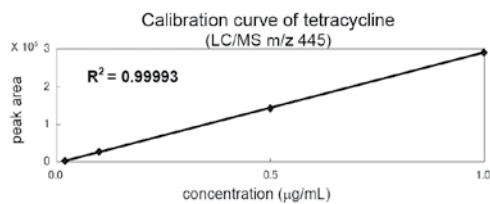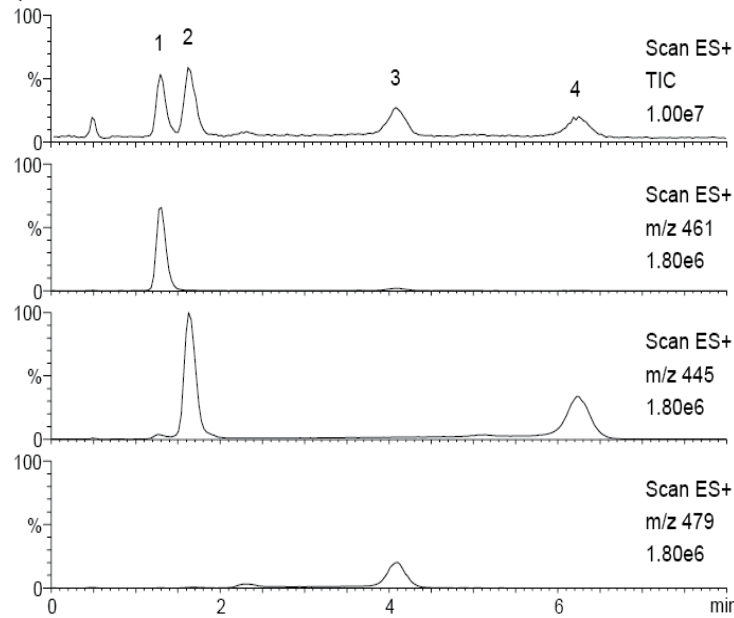


Column: YMC-Triart C18 (5 μm , 12 nm) 50 x 2.0 mm ID
Part-No.: TA12S05-0502WT
Eluent: acetonitrile / water / formic acid (15/85/0.1)
Flow rate: 0.4 ml/min
Temperature: 40 °C
Detection: A) UV at 270 nm
B) ESI positive-mode
Injection: 10 μl (1 $\mu\text{g/ml}$)

Application Data by courtesy YMC Co., Ltd.

22 LC/MS

Analysis of 360 pesticides in a single run

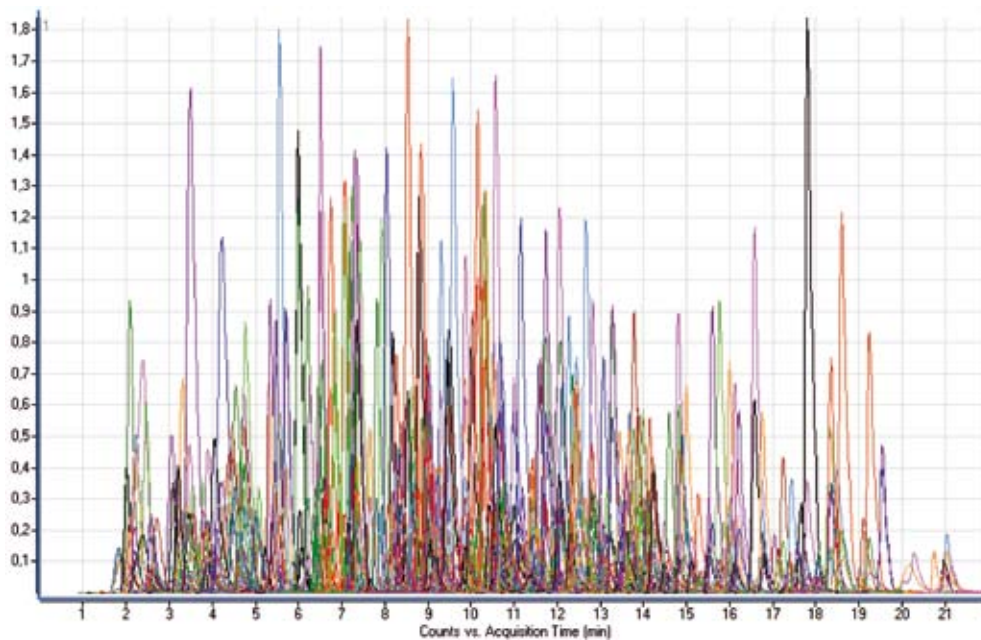


Column: YMC-Triart C18 (3 μ m, 100 x 2.0 mm ID)
Part-No.: TA12503-1002WT
Eluent: A) 5 mM ammonium formate / water
B) 5 mM ammonium formate / methanol
Flow rate: 0.25 ml/min
Temperature: 45 °C

Injection: 5 μ l
Gradient: 0 min: 30% B, 0.1 min: 50% B, 18 min: 100% B,
21 min: 100% B, 21.01 min: 30% B, 29 min: 30% B
Total run time: 30 min
Sample: 100 ng/ml pesticide mix in acetonitrile

by courtesy of: József László
WIREC, WESSLING International Research and Educational Centre Nonprofit Co. (Hungary)



Separation of alkaloids

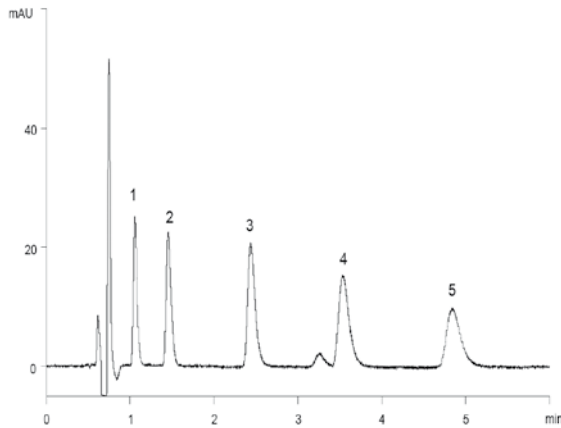

1. Scopolamine
2. Atropine
3. Cinchonine
4. Quinine
5. Dihydroquinine

Column: YMC-Triart C18 (5 μ m, 12 nm)
 50 x 2.0 mm ID
 Part-No.: TA12S05-0502WT
 Eluent: 20 mM CH₃COOH-CH₃COONH₄
 (pH 4.9) / acetonitrile (80/20)
 Flow rate: 0.2 ml/min
 Temperature: 40 °C
 Detection: UV at 220 nm
 Injection: 1 μ l (0.02-0.1 mg/ml)


Barbiturates in human serum

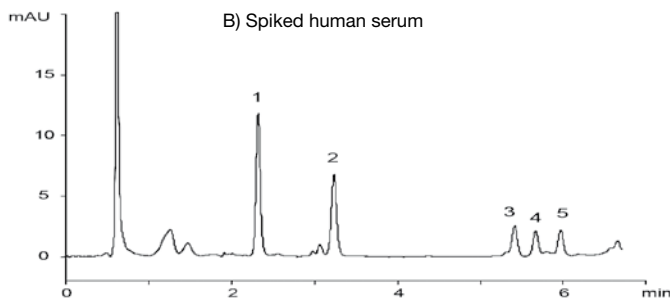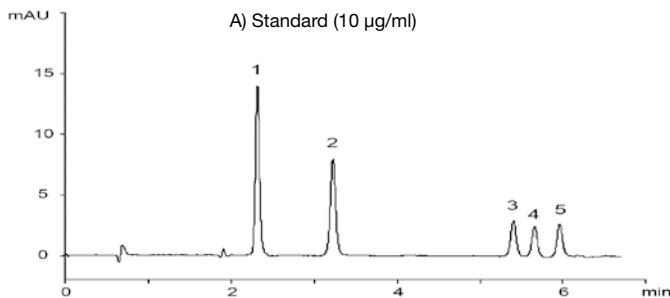

Solid-phase extraction method

YMC Dispo SPE C18 100 mg/1ml

Condition
2 ml methanol 2 ml water
Load
500 μ l spiked human serum solution (each 10 μ g)
Elute
500 μ l methanol/water (85/15)
Dilute
500 μ l 20 mM ammonium formate buffer (pH 9.5)

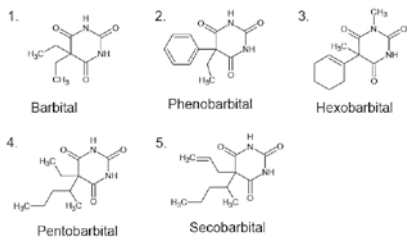


Column: YMC-Triart C18 (5 μ m, 12 nm)
 50 x 2.0 mm ID
 Part-No.: TA12S05-0502WT
 Eluent: A) 20 mM HCOONH₄-NH₃ (pH 9.5)
 B) methanol
 Gradient: 0-90% B (0-7 min)
 Flow rate: 0.2 ml/min
 Temperature: 25 °C
 Detection: UV at 240 nm
 Injection: 1 μ l


Application Data by courtesy YMC Co., Ltd.

Erythromycin at elevated pH and temperature


1. Optimisation of pH

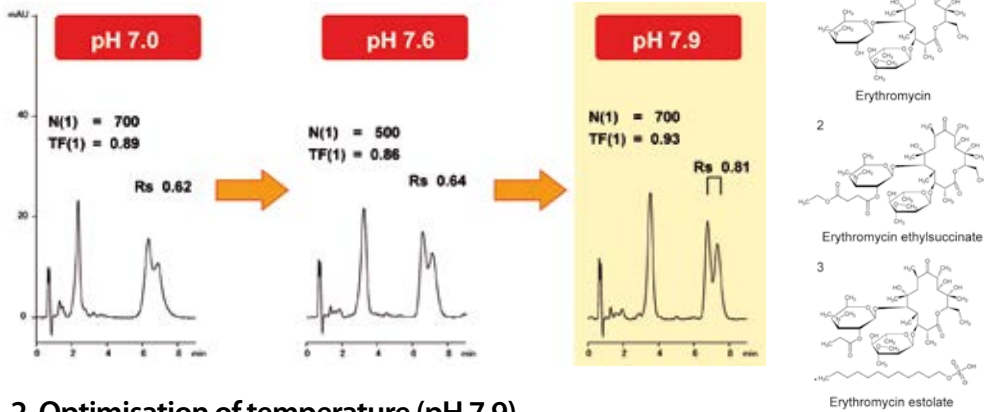





2. Optimisation of temperature (pH 7.9)

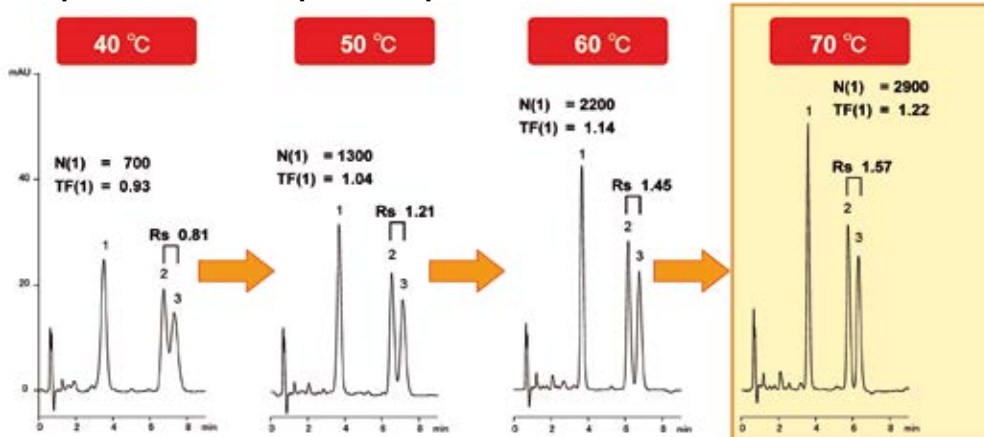

3. Stability test: pH 7.9, 70 °C




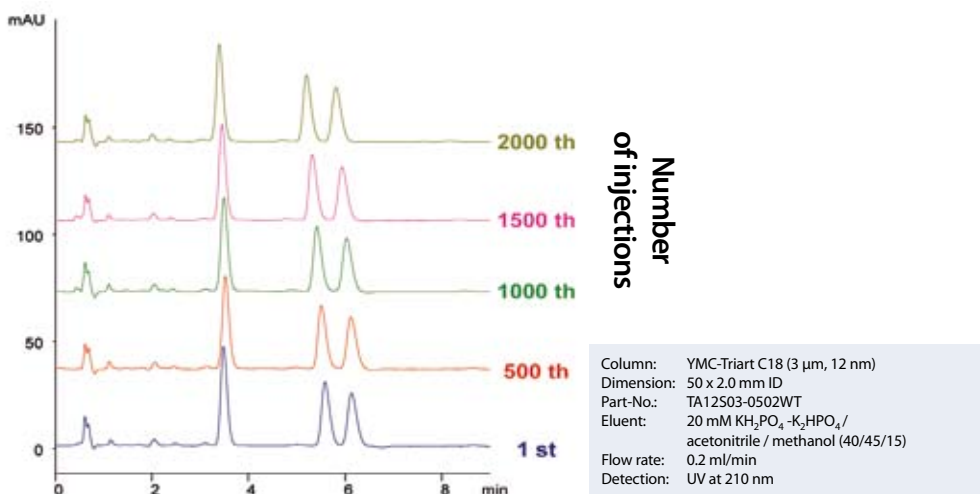


Separation of flavonoids

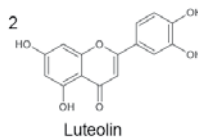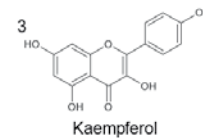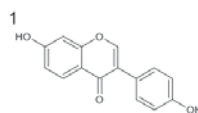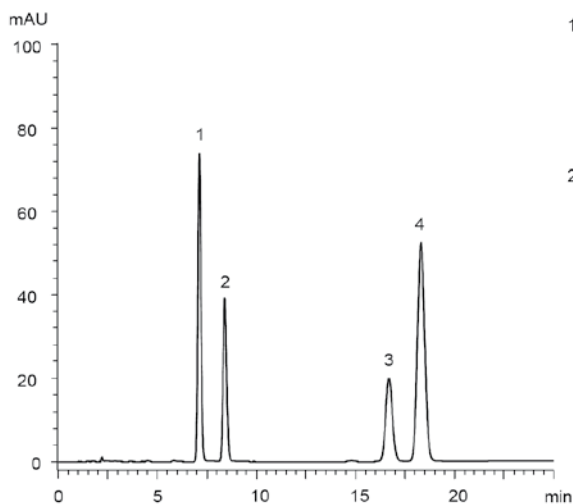

Column: YMC-Triart C18 (5 μ m, 12 nm)
 Dimension: 150 x 3.0 mm ID
 Part-No.: TA12505-1503WT
 Eluent: acetonitrile / 10 mM H₃PO₄ (30/70)
 Flow rate: 0.425 ml/min
 Temperature: 37 °C
 Detection: UV at 280 nm
 Injection: 2 μ l (50 μ g/ml)


Separation of aromatic carboxylic acids

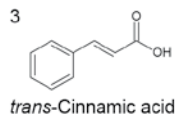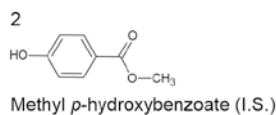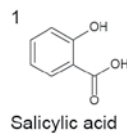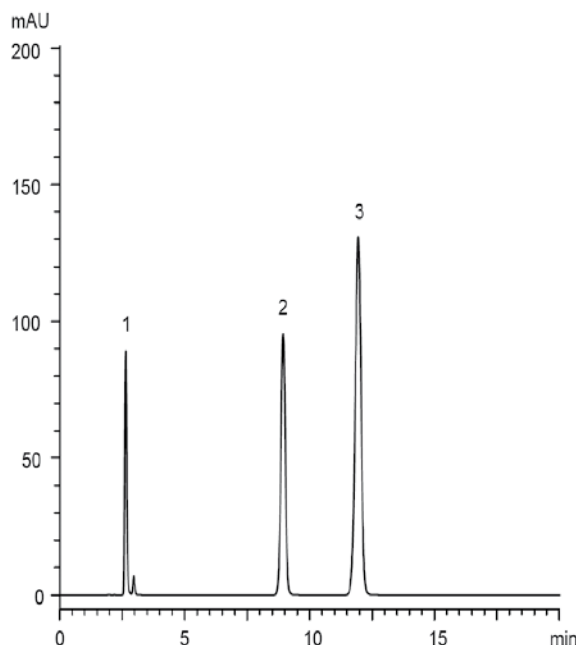


Column: YMC-Triart C18 (5 μ m, 12 nm)
 Dimension: 150 x 3.0 mm ID
 Part-No.: TA12505-1503WT
 Eluent: 10 mM CH₃COOH-CH₃COONH₄ (pH 4.2) / acetonitrile (75/25)
 Flow rate: 0.425 ml/min
 Temperature: 40 °C
 Detection: UV at 254 nm
 Injection: 4 μ l (0.02 ~ 0.3 mg/ml)

Application Data by courtesy YMC Co., Ltd.

Aciclovir syrup and injection

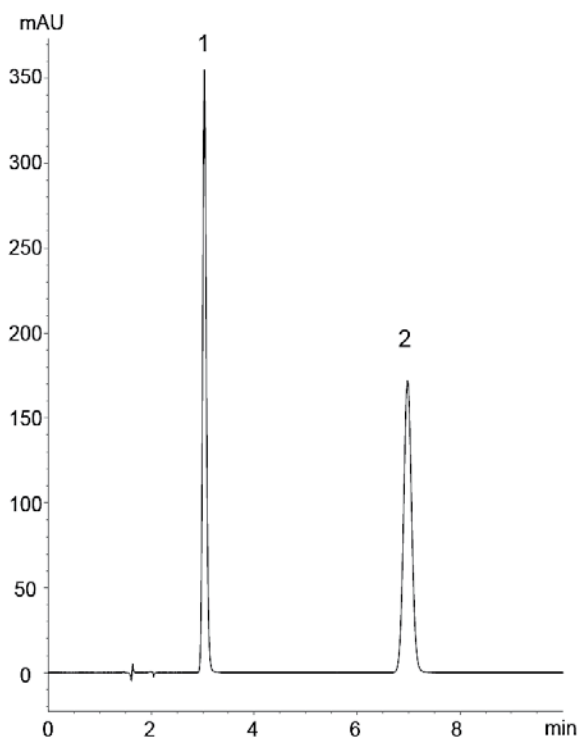


1.

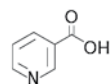

Nicotinic acid



2.

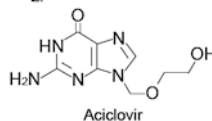


Aciclovir

Column: YMC-Triart C18 (5 μ m, 12 nm)
150 x 4.6 mm ID
Part-No.: TA12S05-1546WT
Eluent: phosphate buffer* / methanol (95/5)
*Dissolve 1.45 g of H₃PO₄ and 25 ml of 1 mol/l CH₃COOH in water to make 900 ml → adjust pH 2.5 by 1 mol/l NaOH → add water to make 1000 ml
Flow rate: 1.0 ml/min
Temperature: 25 °C
Detection: UV at 254 nm
Injection: 20 μ l (0.05 mg/ml, 0.032 mg/ml)


Clindamycin hydrochloride

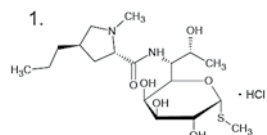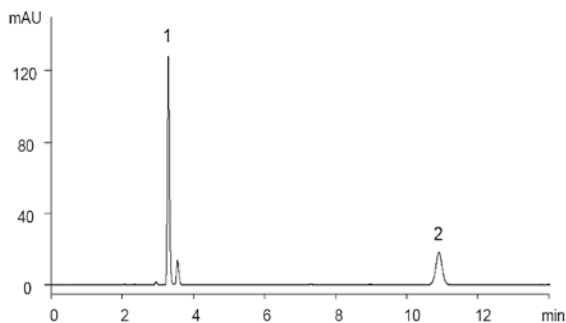

Lincomycin hydrochloride


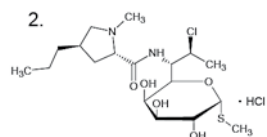


Clindamycin hydrochloride

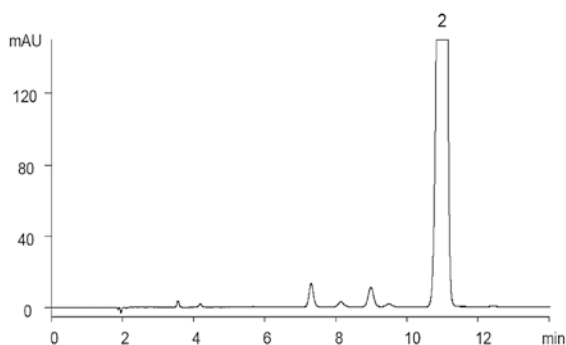


Column: YMC-Triart C18 (5 μ m, 12 nm)
 Dimension: 250 x 4.6 mm ID
 Part-No.: TA12S05-2546WT
 Eluent: 50 mM KH_2PO_4 (pH 7.5 adjusted by 8 M KOH) / acetonitrile (55/45)
 Flow rate: 1.0 ml/min
 Temperature: 25 $^\circ\text{C}$
 Detection: UV at 210 nm
 Injection: 10 μ l

Analysis of amlodipine besilate

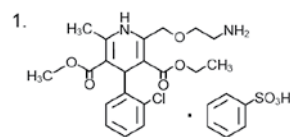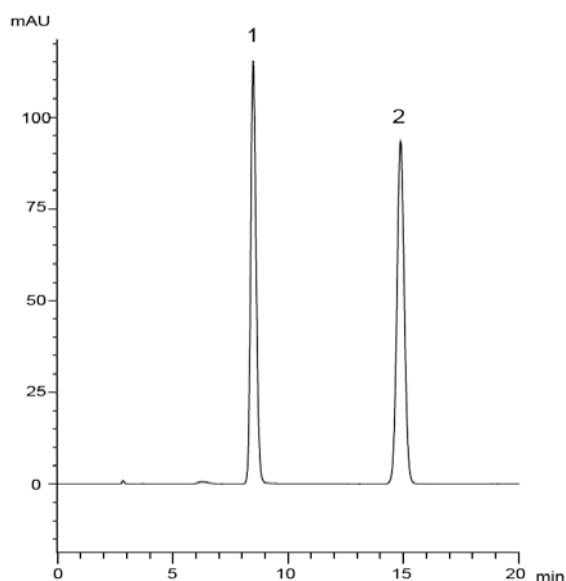

Amlodipine besilate

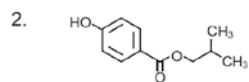


Isobutyl *p*-hydroxybenzoate

Column: YMC-Triart C18 (5 μ m, 12 nm)
 Dimension: 150 x 3.0 mm ID
 Part-No.: TA12S05-1503WT
 Eluent: 10 mM $\text{CH}_3\text{COOH}-\text{CH}_3\text{COONH}_4$ (pH 4.2) / acetonitrile (75/25)
 Flow rate: 0.425 ml/min
 Temperature: 40 $^\circ\text{C}$
 Detection: UV at 254 nm
 Injection: 4 μ l (0.02 ~ 0.3 mg/ml)

Application Data by courtesy YMC Co., Ltd.

Basic drugs

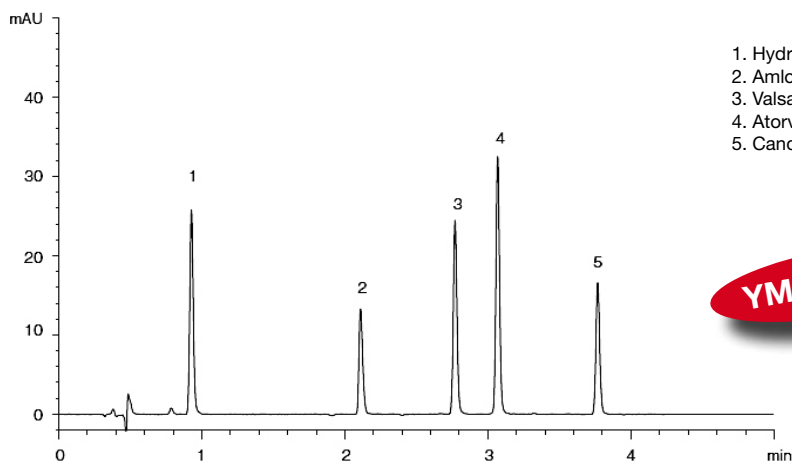



1. Hydrochlorothiazide
2. Amlodipine besilate
3. Valsartan
4. Atorvastatin calcium hydrate
5. Candesartan cilexetil

YMC-Triart C8

Column: YMC-Triart C8 (3 μ m, 12 nm), 50 x 2.0 mm ID
 Part-No.: TO12503-0502WT
 Eluent: A) water / formic acid (100/0.1)
 B) acetonitrile / formic acid (100/0.1)
 10-90% B (0-5 min), 90% B (5-7 min)

Flow rate: 0.4 ml/min
 Temperature: 30 °C
 Detection: UV at 254 nm
 Injection: 2 μ l (10-20 μ g/ml)

Sequential analysis of Carvedilol

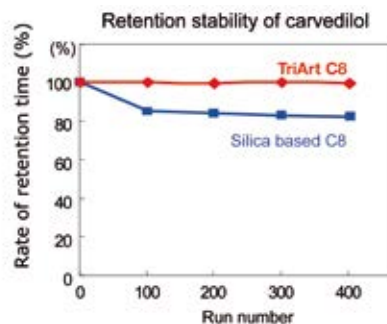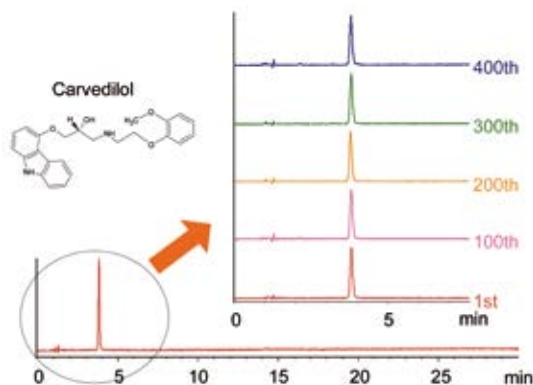


YMC-Triart C8


Column: YMC-Triart C8 (5 μ m, 150 x 2.0 mm ID)
 Part-No.: TO12505-1502WT
 Eluent: phosphate buffer (pH 2.0)* / acetonitrile (65/35)
 *Dissolve 2.72 g of KH_2PO_4 in 900 ml water, adjust pH 2.0 with H_3PO_4 and add water to make 1000 ml
 Flow rate: 0.28 ml/min (adjust the flow rate so that the retention time of carvedilol is about 4 min)
 Temperature: 55 °C
 Detection: UV at 240 nm

No change in retention time is observed even under a high pH and at a elevated temperature.


Soy isoflavones in supplement

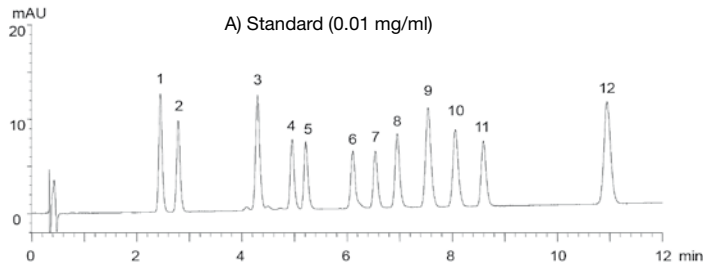

1. Daidzin
2. Glycitin
3. Genistin
4. 6"-O-Malonyldaidzin
5. 6"-O-Malonylglycitin
6. 6"-O-Acetyldaidzin
7. 6"-O-Acetylglycitin
8. 6"-O-Malonylgenistin
9. Daidzein
10. Glycitein
11. 6"-O-Acetylgenistin
12. Genistein

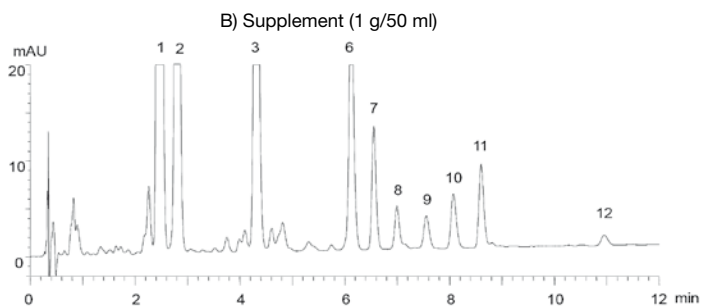

Sample preparation method

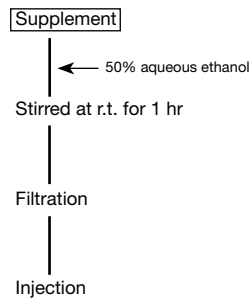

Column: YMC-Triart C18 (3 µm, 12 nm)

50 x 2.0 mm ID

Part-No.: TA12503-0502WT

Eluent: A) acetonitrile / water / HCOOH (10/90/0.1)

B) acetonitrile / water / HCOOH (60/40/0.1)

Gradient: 5-40% B (0-12 min)


Flow rate: 0.4 ml/min

Temperature: 25 °C

Detection: UV at 254 nm

Injection: 2 µl

Separation of water-soluble vitamins

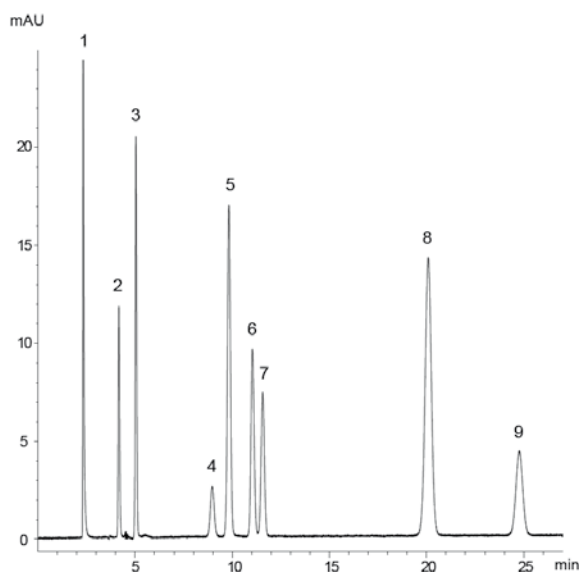

1. Thiamine HCl (Vitamin B₁)
2. Pyridoxine HCl (Vitamin B₆)
3. Nicotinamide
4. Cyanocobalamin (Vitamin B₁₂)
5. L-Ascorbic acid 2-glucoside
6. L-Ascorbic acid (Vitamin C)
7. Erythorbic acid
8. Riboflavin (Vitamin B₂)
9. Nicotinic acid

Column: YMC-Triart C18 (5 µm, 12 nm)

250 x 4.6 mm ID

Part-No.: TA12505-2546WT

Eluent: phosphate buffer* / acetonitrile (90/10)

* Dissolve 1.4 g KH₂PO₄ in 800 ml water

→ add 26 ml 10% TBA:OH

→ adjust pH 5.2 by 20% H₃PO₄

→ add water to make 1000 ml

Flow rate: 0.8 ml/min

Temperature: 40 °C


Detection: UV at 260 nm

Injection: 10 µl (5 µg/ml)

Application Data by courtesy YMC Co., Ltd.

Analysis of anthocyanins and anthocyanidins


Anthocyanins: Indicated in black
 Anthocyanidins: Indicated in red

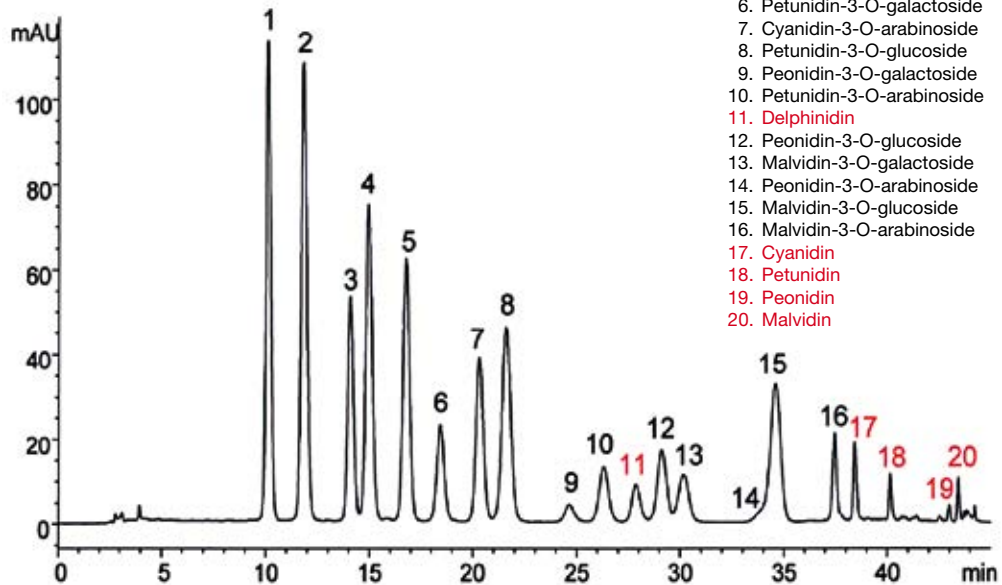

1. Delphinidin-3-O-galactoside
2. Delphinidin-3-O-glucoside
3. Cyanidin-3-O-galactoside
4. Delphinidin-3-O-arabinoside
5. Cyanidin-3-O-glucoside
6. Petunidin-3-O-galactoside
7. Cyanidin-3-O-arabinoside
8. Petunidin-3-O-glucoside
9. Peonidin-3-O-galactoside
10. Petunidin-3-O-arabinoside
11. Delphinidin
12. Peonidin-3-O-glucoside
13. Malvidin-3-O-galactoside
14. Peonidin-3-O-arabinoside
15. Malvidin-3-O-glucoside
16. Malvidin-3-O-arabinoside
17. Cyanidin
18. Petunidin
19. Peonidin
20. Malvidin

Column: YMC-Triart C18 (5 μ m, 12 nm) 250 x 4.6 mm ID
 Part-No.: TA12505-2546WT
 Eluent: A) water / formic acid (90/10)
 B) acetonitrile / methanol / water / formic acid (22.5/22.5/40/10)
 Gradient: 20-28% B (0-30 min), 28-70% B (30-40 min), 100% B (40-45 min)
 Flow rate: 1.0 ml/min
 Temperature: 25 $^{\circ}$ C
 Detection: UV/VIS at 535 nm
 Sample: commercial bilberry powder (1.25 mg/ml)

Peptides (MW 556 - 3,465)

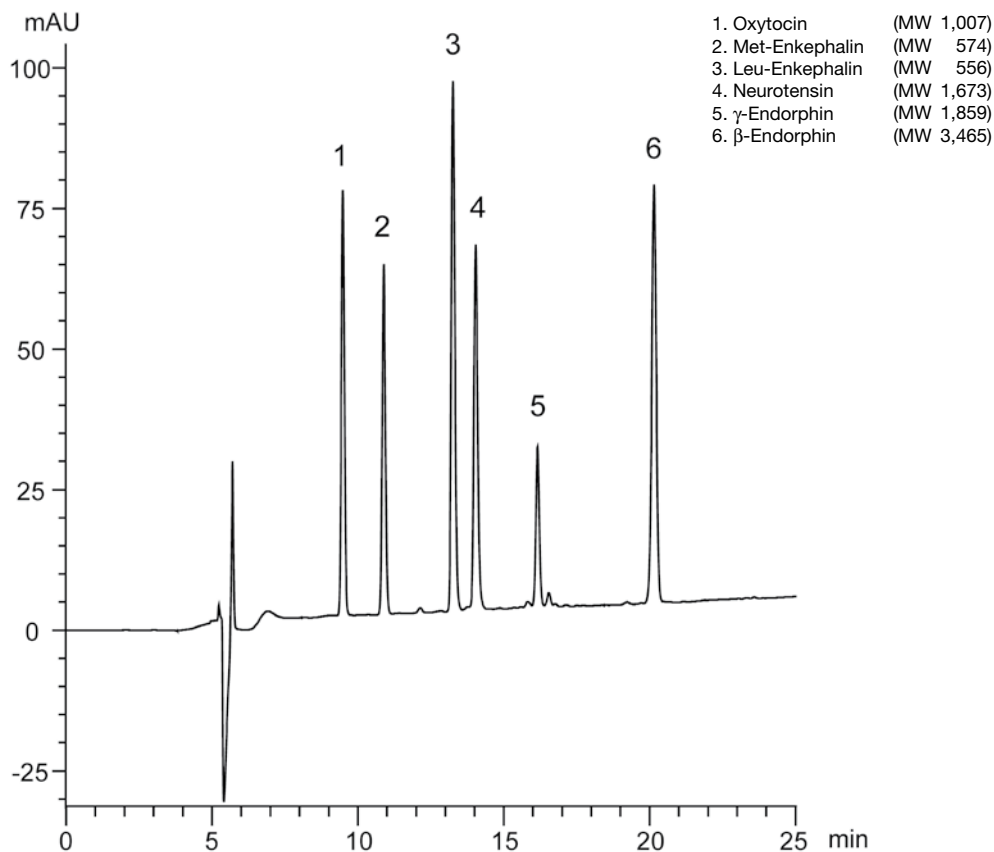

Column: YMC-Triart C18 (5 μ m, 12 nm)
Dimension: 150 x 2.0 mm ID
Part-No.: TA12S05-1502WT
Eluent: A) water / TFA (100/0.1)
B) acetonitrile / TFA (100/0.1)
20-45% B (0-25 min)


Flow rate: 0.2 ml/min
Temperature: 37 $^{\circ}$ C
Detection: UV at 220 nm
Injection: 2 ml (0.075 - 0.25 mg/ml)

Application Data by courtesy YMC Co., Ltd.

32 HILIC

Great durability and reproducibility at high pH

Durability in high pH (pH 11, 50 °C)

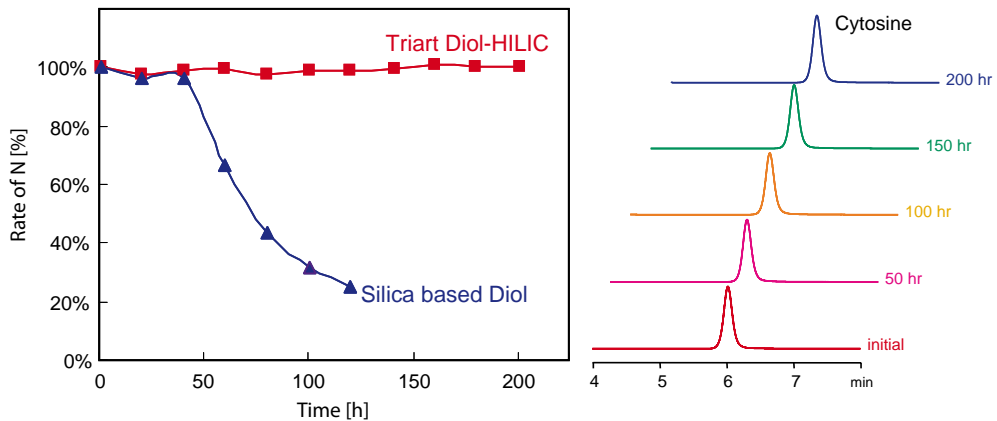


Column: 5 μ m, 150 x 4.6 mm ID
 Part-No.: TDH12S05-1546WT
 Eluent: acetonitrile / water / NH₃ (90/10/0.1) pH 11.3

Flow rate: 1.0 ml/min
 Temperature: 50 °C
 Sample: Cytosine


Triart Diol-HILC offers highly reproducible separations even at high pH and high temperature. The lifetime of Triart Diol-HILIC is much longer than that of conventional silica-based Diol columns.

Nucleosides and bases

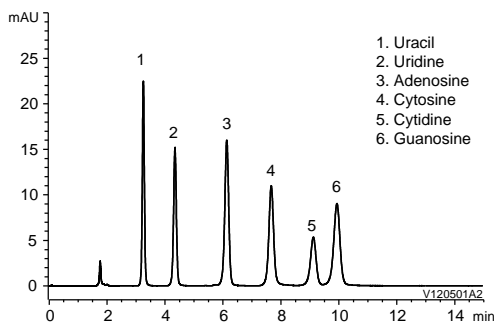


Column: YMC-Triart Diol-HILIC (5 μ m, 12 nm)
 150 x 3.0 mm ID
 Part-No.: TDH12S05-1503WT
 Eluent: 100 mM CH₃COONH₄ / acetonitrile (10/90)
 Flow rate: 0.425 ml/min
 Temperature: 30 °C
 Detection: UV at 254 nm
 Injection: 2 μ l (5 ~ 10 μ g/ml)

Diabetes drugs

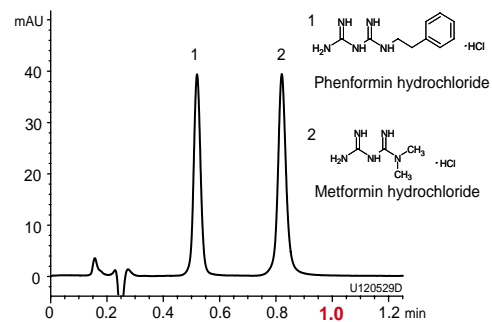

Column: YMC-Triart Diol-HILIC (1.9 μ m, 12 nm)
 50 x 2.0 mm ID
 Part-No.: TDH12SP9-0502PT
 Eluent: 100 mM HCOOH-HCOONH₄ (pH 3.7) / acetonitrile (10/90)
 Flow rate: 0.8 ml/min
 Temperature: 25 °C
 Detection: UV at 235 nm
 Injection: 2 μ l (10 μ g/ml)

Polar and hydrophilic compounds

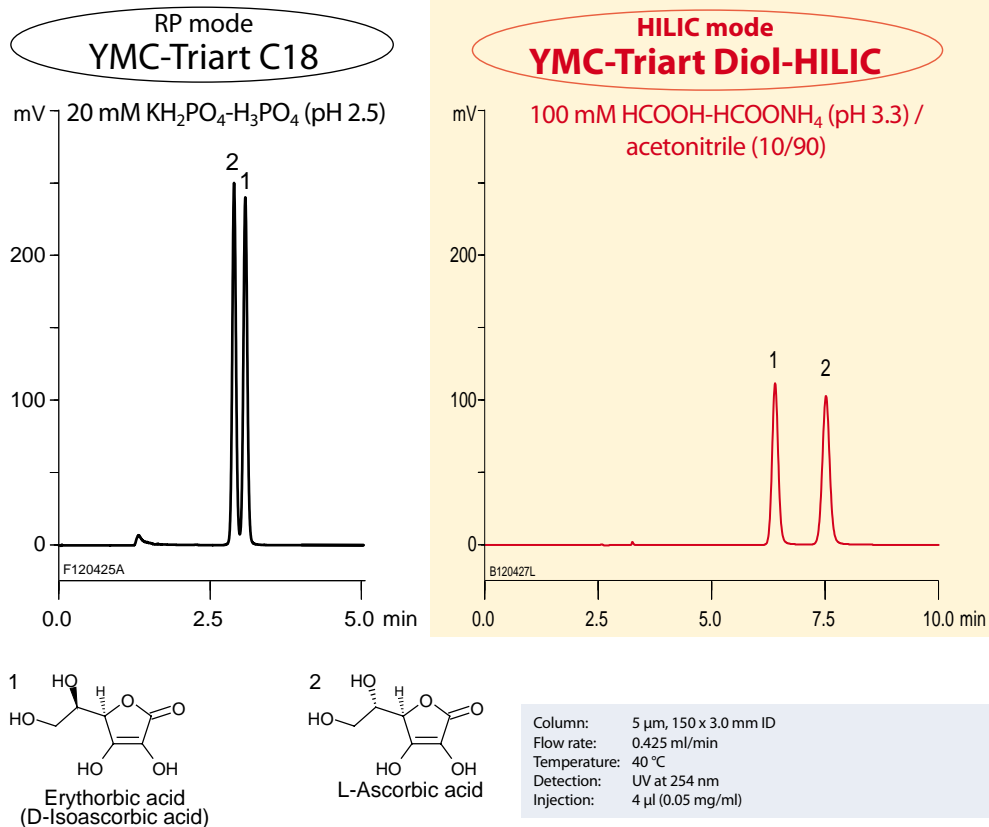


Triart C18 (RP) shows very weak retention and poor resolution of L-ascorbic acid and its stereoisomer (erythorbic acid) even if 100% aqueous mobile phase is used. However, Triart Diol-HILIC shows strong retention and good resolution of these compounds with mobile phase containing 90% organic solvent.

Application Data by courtesy YMC Co., Ltd.


YMC-Triart: Improved quality of particles

Uniform spherical particles

YMC-Triart

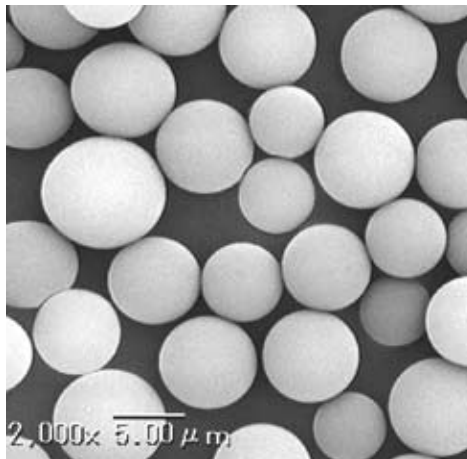


X-Bridge HILIC

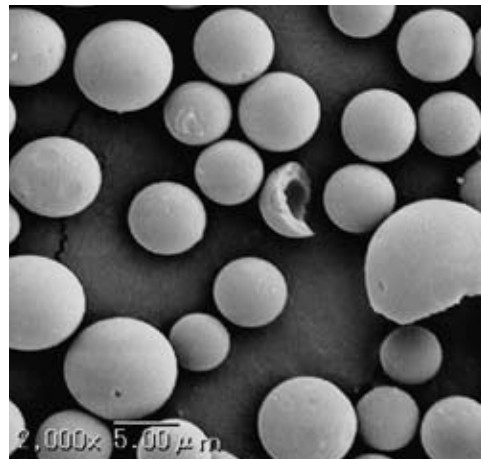

The new uniform spherical particle support is used for YMC-Triart C18 and C8. The particle is produced using **micro-reactor** technology for the granulation process. This results in reduction of the back-pressure and leads to more reproducibility in surface modification.

Low column back-pressure

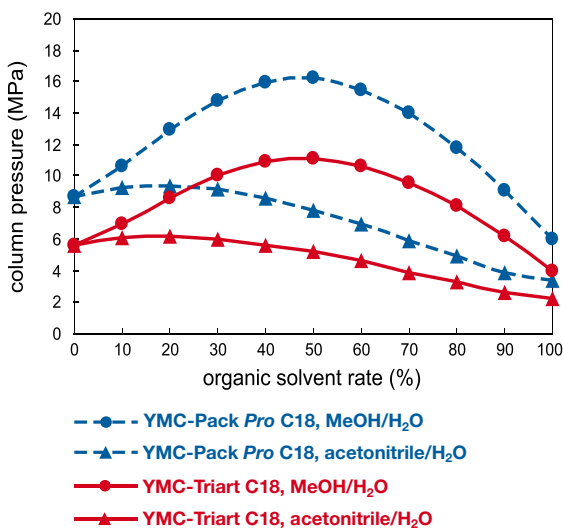

Column: YMC-Triart C18, 5 μm, 150 x 4.6 mm ID
 Part-No.: TA12S05-1546WT
 Eluent: acetonitrile / water or methanol / water
 Flow rate: 1.0 ml/min
 Temperature: 25 °C


The revolutionary production technique, adapted from micro-reactor flow technology, produces a multi-layered silica/organic hybrid stationary phase, with outstanding narrow pore size and particle size distributions which result in low back pressures.

YMC-Triart is designed for use under a wide range of conditions. Elution with higher viscosity methanol (compared with acetonitrile), YMC-Triart generates lower pressure (approx 30% lower than with conventional phases).

Batch-to-Batch reproducibility


Excellent reproducibility of YMC-Triart phases is available even for the analysis for basic and coordination compounds which normally exhibit tailing and adsorption effects.

Basic compounds

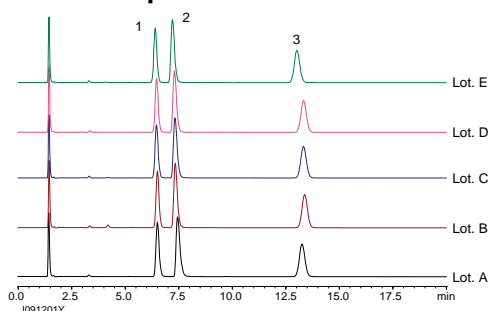


Column: YMC-Triart C18, 5 μ m, 150 x 3.0 mm ID
Part-No.: TA12505-1503WT
Eluent: 20 mM KH_2PO_4 (pH 6.9) / acetonitrile (65/35)
Flow rate: 0.425 ml/min
Temperature: 40 °C
Detection: UV at 235 nm

Coordination compounds

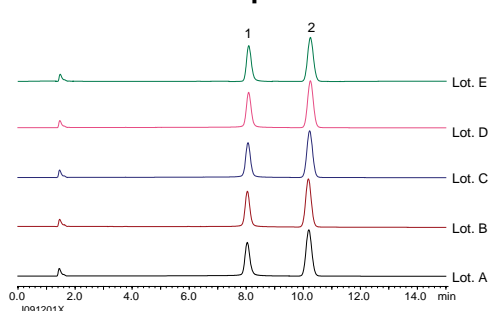

Column: YMC-Triart C18, 5 μ m, 150 x 3.0 mm ID
Part-No.: TA12505-1503WT
Eluent: acetonitrile / 0.1% H_3PO_4 (40/60)
Flow rate: 0.425 ml/min
Temperature: 40 °C
Detection: UV at 254 nm

The reproducibility of packed columns is shown below in terms of theoretical plate number (N) and tailing factor (Tf). YMC-Triart packed columns exhibit a very narrow range of variation.

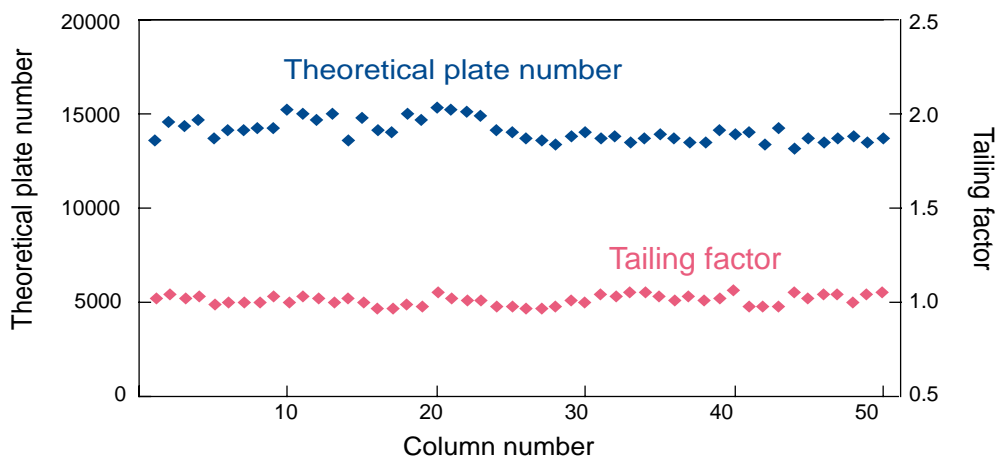


Column: YMC-Triart C18, 5 μ m, 150 x 4.6 mm ID
Part-No.: TA12505-1546WT
Eluent: acetonitrile / water (40/60)
Flow rate: 1.0 ml/min
Temperature: ambient
Sample: butyl benzoate

Application Data by courtesy YMC Co., Ltd.


Multi-stage endcapping

After bonding the alkyl chain, there are highly reactive and less reactive silanols on the surface. In traditional bonding processes, these are reacted with a single capping-compound in one step. However, the highly reactive silanols can be hydrolysed easily which contributes to the poor durability. The less reactive silanols are hard to endcap which results in poor resolution due to peak tailing.

YMC-Triart C18, C8 and Phenyl phases use a new innovation in end capping called "multistage end-capping" for its surface modification process. By using a number of compounds with the different reactivities in successive steps, all silanols can be capped to the maximum extent.

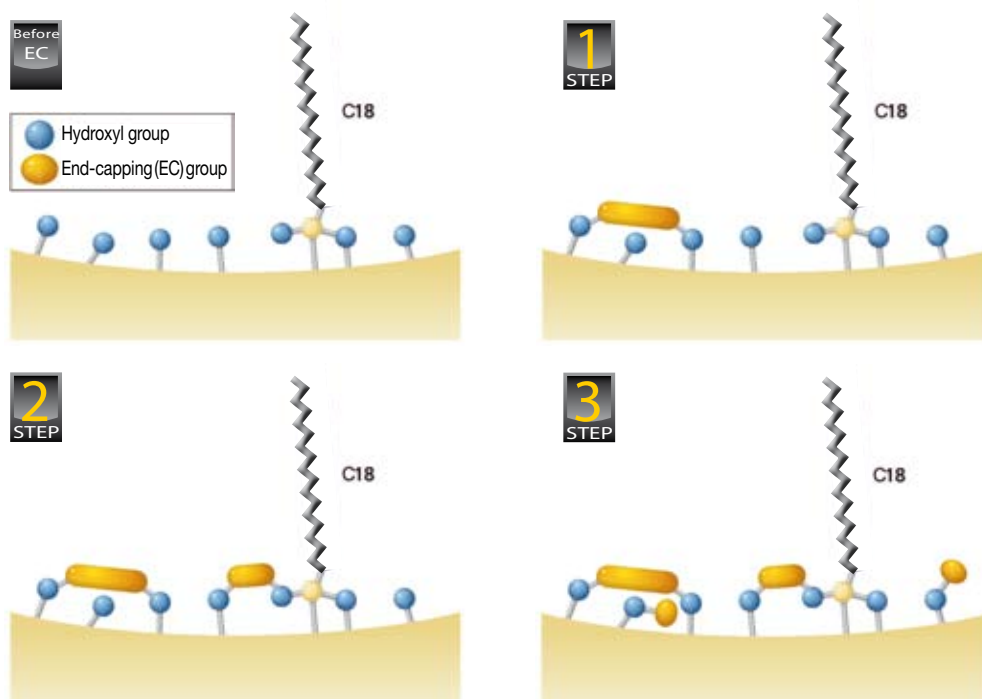


The chromatographic result of a "good" end-capping is demonstrated:

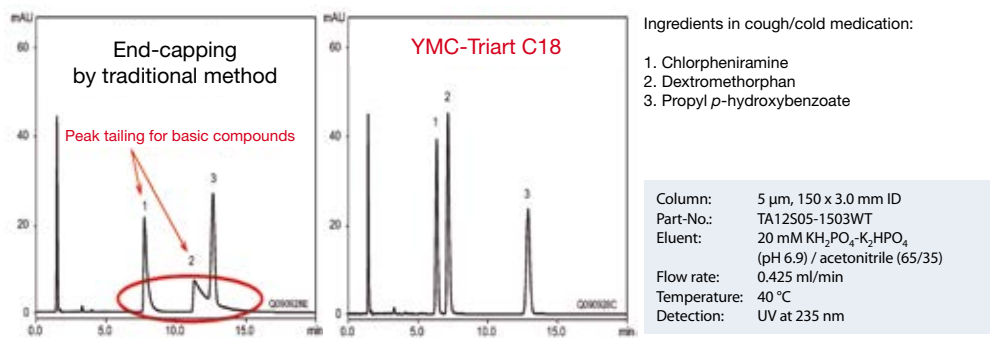

Application Data by courtesy YMC Co., Ltd.

Basic compounds

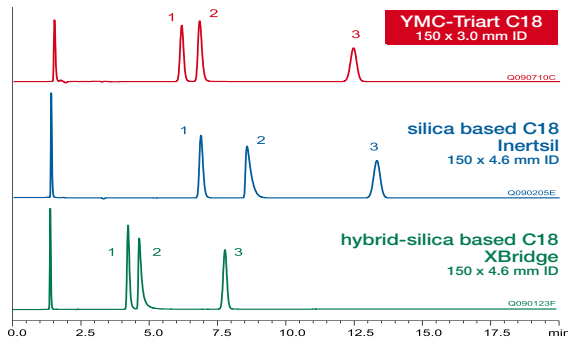


Ingredients in a cough/cold medication

1. Chlorpheniramine
2. Dextromethorphan
3. Propyl p-hydroxybenzoate

Column: 5 μ m, 150 x 3.0 or 150 x 4.6 mm ID
 Eluent: 20 mM KH_2PO_4 - K_2HPO_4 (pH 6.9) / acetonitrile (65/35)
 Flow rate: 0.425 ml/min for 3.0 mm ID
 1.0 ml/min for 4.6 mm ID
 Temperature: 40 °C
 Detection: UV at 235 nm

The innovative surface modification technology results in excellent peak shapes even for basic compounds that often exhibit peak tailing with conventional silica- and hybrid silica-based reversed phase columns.

Acidic compounds

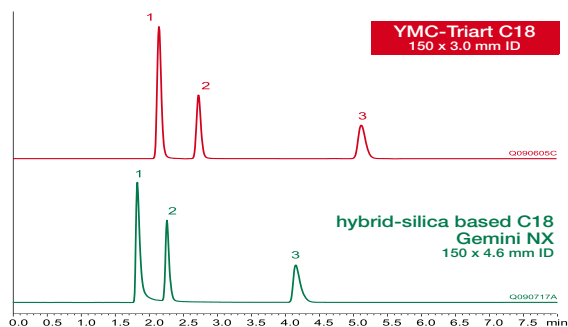


Organic acid

1. Formic acid
2. Acetic acid
3. Propionic acid

Column: 5 μ m, 150 x 3.0 or 150 x 4.6 mm ID
 Eluent: methanol / 0.1% H_3PO_4 (5/95)
 Flow rate: 0.425 ml/min for 3.0 mm ID
 1.0 ml/min for 4.6 mm ID
 Temperature: 37 °C
 Detection: UV at 210 nm

YMC-Triart phases is synthesised using methodology adapted from micro-reactor technology. This technique ensures a reduction of impurities that contribute to peak tailing during the analysis of some types acidic compounds.

Coordinating compounds

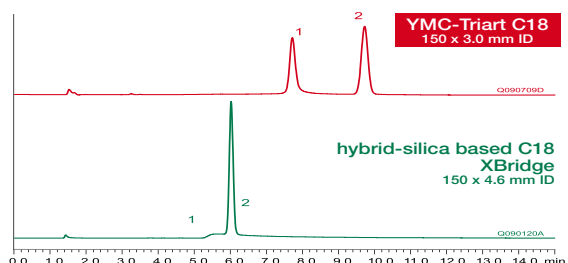



Hinokitiol

1. Hinokitiol
2. Methyl benzoate

Column: 5 μ m, 150 x 3.0 or 150 x 4.6 mm ID
 Eluent: acetonitrile / 0.1% H_3PO_4 (40/60)
 Flow rate: 0.425 ml/min for 3.0 mm ID
 1.0 ml/min for 4.6 mm ID
 Temperature: 40 °C
 Detection: UV at 254 nm

YMC-Triart phases have an extremely low level of metal impurities, much lower than conventional products, ensuring excellent peak shape for coordination compounds.

Comparison of clemastine analysis

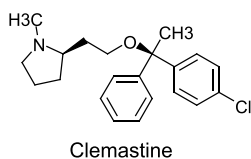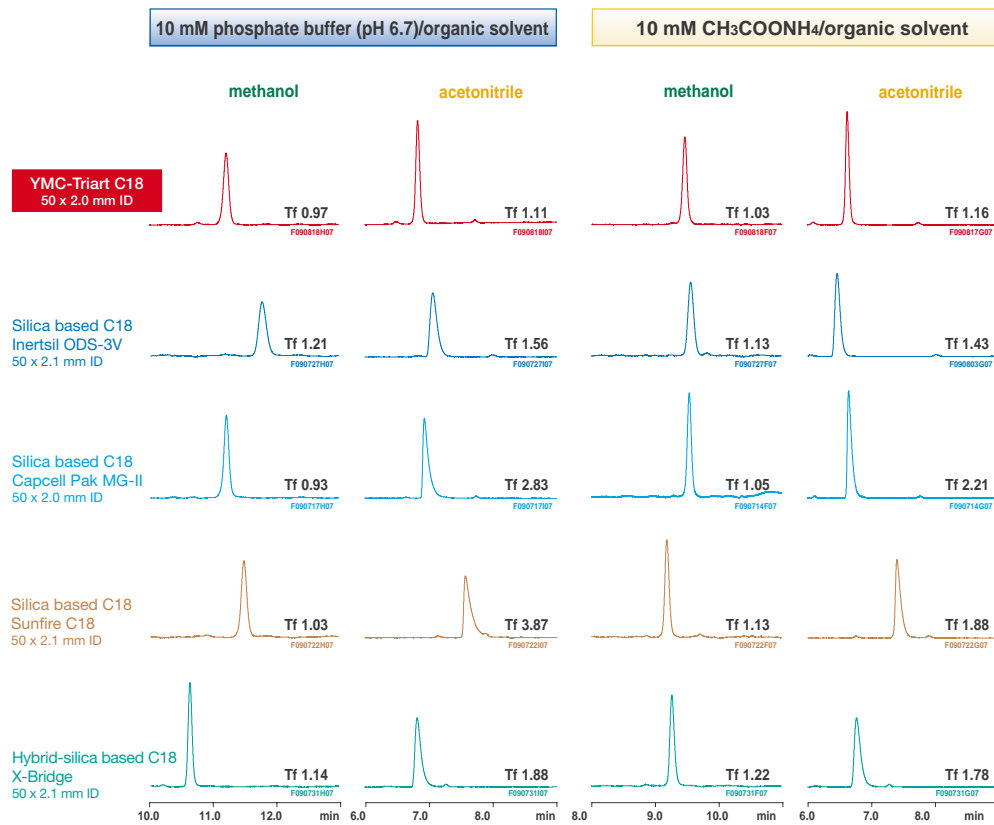

Column: 5 μ m, 50 x 2.0 or 50 x 2.1 mm ID
 Eluent: A) 10 mM KH₂PO₄-K₂HPO₄ (pH 6.7) or 10 mM CH₃COONH₄
 B) methanol or acetonitrile
 5-90% B (0-10 min), 90% B (10-15 min)
 Flow rate: 0.2 ml/min
 Temperature: 25 °C
 Detection: UV at 230 nm

Clemastine is a well-known basic compound which readily exhibits peak tailing with conventional ODS columns. YMC-Triart C18 provides sharp separations with many different buffer/solvent compositions.

Application Data by courtesy YMC Co., Ltd.

Substance index

A	page	C	page	G	page
Acesulfame K	9	4-Chloroacetanilide	18	Genistein	29
4-Acetamidoacetophenone	18	Chlortetracycline	11,21	Genistin	29
2-Acetamidophenol	18,19	trans-Cinnamic acid	25	Glycitein	29
Acetaminophen	18,19	Cinchonine	23	Glycitin	29
Acetanilide	18,19	Clemastine	39	Glycyrrhizin acid ammonium salt	17
Acetic acid	38	Clinamycin HCl	27	Guaiacol	18
6"-O-Acetylaidzin	29	Coordinating compounds	38	Guanosine	32
6"-O-Acetylgenistin	29	Copper 8-quinolinolate	7	H	
6"-O-Acetylglycitin	29	Cyanidin	30	Halosulfuronmethyl	7
Acetylsalicylic acid	19	Cyanidin-3-O-arabinoside	30	Herbicides	7
Aciclovir	26	Cyanidin-3-O-galactoside	30	Hexobarbital	23
Adenosine	32	Cyanidin-3-O-glucoside	30	Hinokitiol	38
Adrenaline hydrochloride (A)	15	Cyanocobalamin	29	Human insulin	10
Agrochemicals	7	Cyclamate Na	9	Homovanillic acid (HVA)	15
Acidic compounds	38	Cytidine	32	Hydrochlorothiazide	28
Alkaloids	23	Cytosine	32	2-Hydroxyacetophenone	18
4-Aminophenone	18	D		5-Hydroxyindoleacetic acid (5HIAA)	15
Amitriptyline	16	Daidzein	25,29	5-Hydroxytryptamine hydrochloride (5HT)	15
Amlodipine besilate	27,28	Daidzin	29	I	
Angiotensin I	10	Delphinidin	30	Imipramine	16
Angiotensin II	10	Delphinidin-3-O-arabinoside	30	Insulin	10
Angiotensin III	10	Delphinidin-3-O-galactoside	30	D-Isoascorbic acid	33
Anthocyanidins	30	Delphinidin-3-O-glucoside	30	Isobutyl p-hydroxybenzoate	27
Anthocyanins	30	Dextromethorphan	6,8,37,38	K	
Antidepressants	16	Dihydroquinine	23	Kaempferol	25
Aromatic carboxylic acids	25	3,4-Dihydroxymandelic acid (DOMA)	15	L	
Artificial sweeteners	9	3,4-Dihydroxyphenylacetic acid (DOPAC)	15	Lincomycin HCl	27
L-Ascorbic acid	29,33	3,4-Dihydroxyphenylalanine (DOPA)	15	Luteolin	25
L-Ascorbic acid 2-glucoside	29	1,2-Dimethoxyl benzene	18	M	
Asulam	7	Diquat	11	Macrolide antibiotics	10
Atenolol	7	Dopamine hydrochloride (DA)	15	Maleic acid sodium salt	17
Atorvastatin calcium hydrate	28	Doxycycline	21	6"-O-Malonyldaidzin	29
Atropine	23	E		6"-O-Malonylgenistin	29
Azoxystrobin	7	Endorphins	31	6"-O-Malonylglycitin	29
B		Enkephalins	31	Malvidin	30
Baicalein	25	Epinephrine hydrochloride	15	Malvidin-3-O-arabinoside	30
Barbital	23	Erythorbic acid	29,33	Malvidin-3-O-galactoside	30
Barbiturates	23	Erythromycin	10,24	Malvidin-3-O-glucoside	30
Benzethonium chloride	17	Erythromycin estolate	24	Mecoprop	7
Benzodiazepine derivatives	20	Erythromycin ethylsuccinate	24	Metformin hydrochloride	32
Basic compounds	38	Etizolam	20	3-Methoxy-4-hydroxyphenylglycol (MHPG)	15
Basic drugs	28	F		3-Methoxytyramine hydrochloride (3MT)	15
Betablockers	11	Fexofenadine hydrochloride	12,13	Methyl benzoate	18,38
Bovine insulin	10	Flavonoids	25	4-Methylcatechol	16
Brilliant blue FCF	14	Flazasulfuron	7	Methyl p-hydroxybenzoate	25
C		Formic acid	38		
Caffeine	19	Fungicides	7		
Candesartan cilexetil	28				
Carvedilol	28				
Catecholamines	15				
Chlorphenamine	8				
Chlorpheniramine	6,17,37,38				

Substance index

M	page	P	page	S	page
2-Methylresorcinol	16	Petunidin	30	Stevioside hydrate	17
2,5-Methylresorcinol	16	Petunidin-3-O-arabinoside	30	Sulpha drugs	10
Metoprolol	9	Petunidin-3-O-galactoside	30	Sulphamerazin	10
		Petunidin-3-O-glucoside	30	Sulphamethoxazole	10
N		Phenacetine	19	Sulphathiazole	10
Nadolol	11	Phenobarbital	23		
Naphazolin HCl	17	Phenoformin hydrochloride	32	T	
Neostigmine methylsulfate	17	Phenol	18	Tetracycline antibiotics	11,21
Neurotensin	31	Pindolol	11	Tetrahydrozoline HCl	17
Nicotinamide	29	Procine insulin	10	Thiamine HCl	29
Nicotinic acid	26,29	Propranolol	11,36	Thiram	7
4-Nitrocatechol	16	Propionic acid	38	Triazolam	20
4-Nitrophenol	18	Propyl paraben	6,8	Triclopyr	7
Noradrenaline hydrochloride (NA)	15	Propyl p-hydroxybenzoate	37,38	Tryptophan (Trp)	15
Norephedrine bitartate	16	Pyridoxine HCl	17,29	Tyrosine (Tyr)	15
Norepinephrine hydrochloride	15	Pyrocatechol	16,18		
Nortriptyline	16			U	
O		Q		Uracil	11,32
Organic acids	38	Quinine	32	Uridine	32
Oxine-copper	7	R		V	
Oxytetracycline	11,21	Rebaudioside A	17	Valsartan	28
Oxytocin	31	Resorcinole	16	Vanillylmandelic acid (VMA)	15
		Riboflavin	29	Verapamil	36
P		S		Vitamin B1	29
Paracetamol	18	Saccharin	9	Vitamin B2	29
Paraquat	11	Salicylic acid	19,25	Vitamin B6	29
Pentobarbital	23	Scopolamine	23	Vitamin B12	29
Peonidin	30	Secobarbital	23	Vitamin C	29
Peonidin-3-O-arabinoside	30	Serotonin hydrochloride	15	W	
Peonidin-3-O-galactoside	30	Siduron	7	Water-soluble vitamins	29,32
Peonidin-3-O-glucoside	30	Soy isoflavones	29		
Peptides	31	Spiramycin	10		
Pesticides	7,22				

The following brands, trademarks, or service marks are the property of the listed company and/or its subsidiaries. Every effort has been taken to ensure this list is accurate at the time of printing this brochure.

SureFit is trademark of IDEX Corporation.

Capcell Pak MG-II, Capcell Pak C18 AQ are trademarks of Shiseido Co., Ltd.

Intertsil ODS-3V, Intertsil ODS-3, Intertsil ODS-4 are trademarks of GL Sciences, inc..

Acquity UPLC BEH C18, Acquity UPLC CSH C18, Atlantis dC18, Atlantis T3, SunFire C18, Symmetry C18, X-Bridge C18, XTerra MS C18 are trademarks of Waters Corp.

Gemini, Gemini NX, Kinetex C18 are trademarks of Phenomenex, Inc.

Hypersil GOLD is a trademark of Thermo Fisher Scientific.

Kromasil Eternity is a trademark of EKA Chemicals AB.

TSK-gel ODS-100S is a trademark of Tosoh Bioscience LLC.

Zorbax Extend-C18, Zorbax StableBond, Zorbax Eclipse Plus are trademarks of Agilent Technologies, Inc..

L-column is a trademark of Chemicals Evaluation and Research Institute.

Ordering information

YMC-Triart 1.9 µm UHPLC columns

Phase	Column ID (mm)	Column length (mm)					Guard cartridges* with 5 mm length (pack of 3)
		30	50	75	100	150	
C18	2.0	TA12SP9-0302PT	TA12SP9-0502PT	TA12SP9-L502PT	TA12SP9-1002PT	TA12SP9-1502PT	TA12SP9-E5Q1CC
	3.0	—	TA12SP9-0503PT	TA12SP9-L503PT	TA12SP9-1003PT	TA12SP9-1503PT	TA12SP9-E503CC
C8	2.0	TO12SP9-0302PT	TO12SP9-0502PT	TO12SP9-L502PT	TO12SP9-1002PT	TO12SP9-1502PT	TO12SP9-E5Q1CC
	3.0	—	TO12SP9-0503PT	TO12SP9-L503PT	TO12SP9-1003PT	TO12SP9-1503PT	TO12SP9-E503CC
Phenyl	2.0	—	TPH12SP9-0502PT	—	TPH12SP9-1002PT	—	TPH12SP9-E5Q1CC
	3.0	—	—	—	—	—	TPH12SP9-E503CC
PPF	2.0	—	TPF12SP9-0502PT	—	TPF12SP9-1002PT	—	TPF12SP9-E5Q1CC
	3.0	—	—	—	—	—	TPF12SP9-E503CC
HILIC	2.0	—	TDH12SP9-0502PT	—	TDH12SP9-1002PT	—	TDH12SP9-E5Q1CC
	3.0	—	TDH12SP9-0503PT	—	TDH12SP9-1003PT	—	TDH12SP9-E503CC

*Guard cartridge holder required, part no. XPCHUHP

YMC-Triart 3 µm analytical columns

Phase	Column ID (mm)	Column length (mm)					Guard cartridges* with 10 mm length (pack of 5)
		50	75	100	150	250	
C18	2.0	TA12S03-0502WT	TA12S03-L502WT	TA12S03-1002WT	TA12S03-1502WT	—	TA12S03-01Q1GC
	3.0	TA12S03-0503WT	TA12S03-L503WT	TA12S03-1003WT	TA12S03-1503WT	—	TA12S03-0103GC
	4.6	TA12S03-0546WT	TA12S03-L546WT	TA12S03-1046WT	TA12S03-1546WT	TA12S03-2546WT	TA12S03-0104GC
C8	2.0	TO12S03-0502WT	TO12S03-L502WT	TO12S03-1002WT	TO12S03-1502WT	—	TO12S03-01Q1GC
	3.0	TO12S03-0503WT	TO12S03-L503WT	TO12S03-1003WT	TO12S03-1503WT	—	TO12S03-0103GC
	4.6	TO12S03-0546WT	TO12S03-L546WT	TO12S03-1046WT	TO12S03-1546WT	TO12S03-2546WT	TO12S03-0104GC
Phenyl	2.0	TPH12S03-0502WT	—	—	TPH12S03-1502WT	—	TPH12S03-01Q1GC
	3.0	—	—	TPH12S03-1003WT	TPH12S03-1503WT	—	TPH12S03-0103GC
	4.6	TPH12S03-0546WT	TPH12S03-L546WT	TPH12S03-1046WT	TPH12S03-1546WT	—	TPH12S03-0104GC
PPF	2.0	TPF12S03-0502WT	—	—	TPF12S03-1502WT	—	TPF12S03-01Q1GC
	3.0	—	—	TPF12S03-1003WT	TPF12S03-1503WT	—	TPF12S03-0103GC
	4.6	TPF12S03-0546WT	TPF12S03-L546WT	TPF12S03-1046WT	TPF12S03-1546WT	—	TPF12S03-0104GC
HILIC	2.0	TDH12S03-0502WT	—	TDH12S03-1002WT	TDH12S03-1502WT	—	TDH12S03-01Q1GC
	3.0	TDH12S03-0503WT	—	TDH12S03-1003WT	TDH12S03-1503WT	—	TDH12S03-0103GC
	4.6	—	—	TDH12S03-1046WT	TDH12S03-1546WT	—	TDH12S03-0104GC

*Guard cartridge holder required, part no. XPGCH-Q1

YMC-Triart 5 µm analytical columns

Phase	Column ID (mm)	Column length (mm)					Guard cartridges* with 10 mm length (pack of 5)
		50	75	100	150	250	
C18	2.0	TA12S05-0502WT	TA12S05-L502WT	TA12S05-1002WT	TA12S05-1502WT	—	TA12S05-01Q1GC
	3.0	TA12S05-0503WT	TA12S05-L503WT	TA12S05-1003WT	TA12S05-1503WT	—	TA12S05-0103GC
	4.6	TA12S05-0546WT	TA12S05-L546WT	TA12S05-1046WT	TA12S05-1546WT	TA12S05-2546WT	TA12S05-0104GC
C8	2.0	TO12S05-0502WT	TO12S05-L502WT	TO12S05-1002WT	TO12S05-1502WT	—	TO12S05-01Q1GC
	3.0	TO12S05-0503WT	TO12S05-L503WT	TO12S05-1003WT	TO12S05-1503WT	—	TO12S05-0103GC
	4.6	TO12S05-0546WT	TO12S05-L546WT	TO12S05-1046WT	TO12S05-1546WT	TO12S05-2546WT	TO12S05-0104GC
Phenyl	2.0	—	—	—	TPH12S05-1502WT	—	TPH12S05-01Q1GC
	3.0	TPH12S05-0503WT	—	—	—	—	TPH12S05-0103GC
	4.6	TPH12S05-0546WT	—	TPH12S05-1046WT	TPH12S05-1546WT	TPH12S05-2546WT	TPH12S05-0104GC
PPF	2.0	—	—	—	TPF12S05-1502WT	—	TPF12S05-01Q1GC
	3.0	TPF12S05-0503WT	—	—	—	—	TPF12S05-0103GC
	4.6	TPF12S05-0546WT	—	TPF12S05-1046WT	TPF12S05-1546WT	TPF12S05-2546WT	TPF12S05-0104GC
HILIC	2.0	—	—	—	TDH12S05-1502WT	—	TDH12S05-01Q1GC
	3.0	—	—	—	TDH12S05-1503WT	—	TDH12S05-0103GC
	4.6	—	—	—	TDH12S05-1546WT	TDH12S05-2546WT	TDH12S05-0104GC

*Guard cartridge holder required, part no. XPGCH-Q1

Ordering Information

YMC-Triart, 12 nm, 5 µm in ACTUS high-throughput semipreparative hardware

Phase	Column ID (mm)	Column length (mm)				
		50	75	100	150	250
C18	20.0	TA12S05-0520WX	—	TA12S05-1020WX	TA12S05-1520WX	TA12S05-2520WX
	30.0	TA12S05-0530WX	TA12S05-L530WX	TA12S05-1030WX	TA12S05-1530WX	TA12S05-2530WX
C8	20.0	TO12S05-0520WX	—	TO12S05-1020WX	TO12S05-1520WX	TO12S05-2520WX
	30.0	TO12S05-0530WX	TO12S05-L530WX	TO12S05-1030WX	TO12S05-1530WX	TO12S05-2530WX

YMC-Triart, preparative bulk media

YMC-Triart C18-S			YMC-Triart C8-S		
Pore size (nm)	Particle size (µm)	Product Code	Pore size (nm)	Particle size (µm)	Product Code
12	10	TAS12S11	20	10	TOS20S11
	15	TAS12S16		15	TOS20S16
	20	TAS12S21		20	TOS20S21
	—	—		—	—

Available in pack sizes 100 g, 500 g, 1 kg, 5 kg


Please inquire for the corresponding brochures:

YMC-Triart for AQ applications (16 p.)




YMC-Triart Prep Brochure (12 p.)



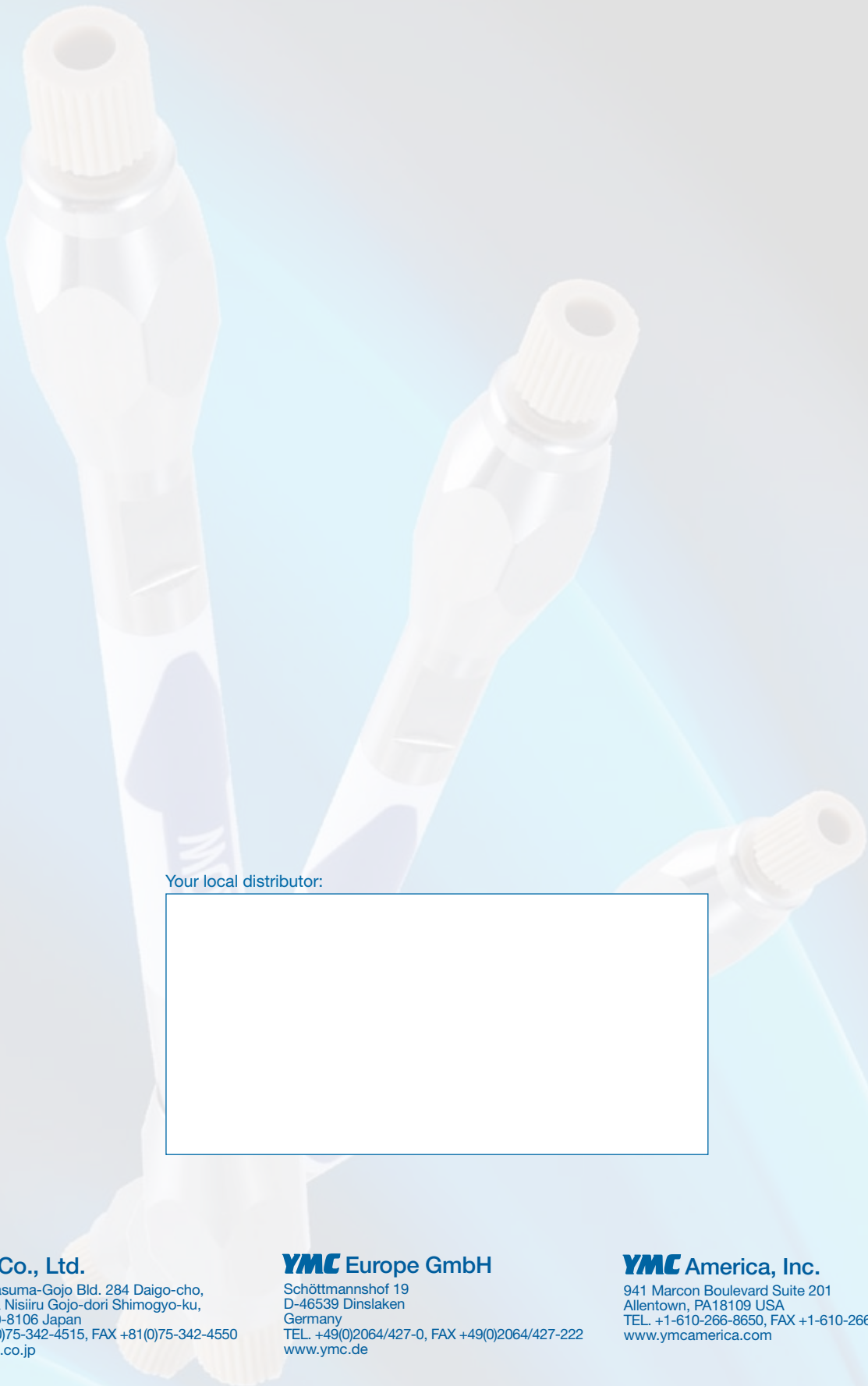

Your local distributor:

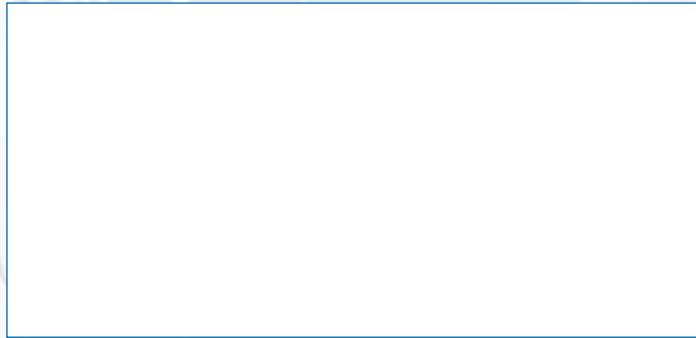

YMC Co., Ltd.

YMC Karasuma-Gojo Bld. 284 Daigo-cho,
Karasuma Nisuru Gojo-dori Shimogyo-ku,
Kyoto 600-8106 Japan
TEL. +81(0)75-342-4515, FAX +81(0)75-342-4550
www.ymc.co.jp

YMC Europe GmbH

Schöttmannshof 19
D-46539 Dinslaken
Germany
TEL. +49(0)2064/427-0, FAX +49(0)2064/427-222
www.ymc.de

YMC America, Inc.

941 Marcon Boulevard Suite 201
Allentown, PA18109 USA
TEL. +1-610-266-8650, FAX +1-610-266-8652
www.ymcamerica.com

YMC Taiwan Co., Ltd.

3F, No. 1353, Zhongzheng Rd.,
Taoyuan City, Taoyuan Country 330,
Taiwan (R.O.C.)
TEL. +886-3-2150-630, FAX +886-3-2150-286
www.ymctaiwan.com

YMC Korea Co., Ltd.

#208, Owners Tower, 16-5, Sunae-dong,
Bundang-gu, Seongnam-si, Gyeonggi-do,
463-825 Korea
TEL. +82-31-716-1631, FAX +82-31-716-1630
www.ymckorea.com

YMC India Pvt. Ltd.

A-154, Eros Boulevard Hilton Hotel Complex Plot
No.13-B, District Centre Mayur Vihar Phase-I
New Delhi-110091 India
TEL. +91-11-45041601;45041701, FAX +91-11-45041901
www.ymcindia.com

YMC Co., Ltd. Shanghai Rep. Office

Far East International Plaza A2404
No. 319 Xianxia Road, Shanghai 200051
P.R. China
TEL. +86-21-6235-1388, FAX +86-21-6235-1398
www.ymcchina.com

YMC Singapore Tradelinks Pte. Ltd.

6 Indus Road, #13-09 Emerald Park, Tower 1,
Singapore 169588
TEL.+65-9023-7617
www.ymc.sg